CUVEIIE BUYING GUIDE

A comprehensive list of premium quality cuvettes for a wide variety of applications and instruments.



- Spectrophotometer cuvettes
- Fluorometer cuvettes
- Cuvette accessories
- Tablet dissolution & flow cells
- Sub-micro volume cuvettes
- Dye laser cells
- Plates & discs
- Microfluidic chips
- Flow channel cells

- · Gas cells
- · Custom manufacturing
- UV/VIS/NIR spectrophotometer calibration standards
- Fluorescence calibration standards
- Evacuable dies
- Refractometers
- and much more...



PRECISION QUARTZ GLASS MANUFACTURING

TABLE OF CONTENTS

ORDERING INFO	NG INFORMATION (OEMs) 3 DRMATION 4 NICAL INFORMATION 5	I	NDEX BY TYPE	76
SECTION 1	SPECTROPHOTOMETER CUVETTES	S	SECTION 7	MICROFLUIDIC CHIPS
1-1 1-2 1-3 1-4	Standard Cuvette Cylindrical Cuvettes Micro, Semi-Micro, and Self Masking Cuvettes Sub-Micro Cuvettes		7-1 7-2	Standard Microfluidic Chips Custom Microfluidic Chips
1-5 1-6	Ultra Micro Cuvette Screw Cap Cuvettes	S	SECTION 8	FLOW CHANNEL CELLS
1-7 1-8 1-9 1-10 1-11	Flow Cells/Tablet Dissolution Cells Micro and Self Masking Flow Cells Water Jacketed Cuvettes Anaerobic Cuvettes Tandem (Divided) Cuvettes		8-1 8-2	Standard Flow Channel Cell Custom Flow Channel Cells
1-12 1-13	Demountable Cuvettes Cuvettes with Quartz to Glass Grading	S	SECTION 9	OTHER ITEMS
1-14 1-15	Cryogenic Cuvettes Tablet Dissolution Cells		9-1 9-2 9-3	Evacuable Dies and Accessories Cuvette Holders for IR and FTIR Instruments Gas Cells
SECTION 2	FLUOROMETER CUVETTES			
2-1 2-2 2-3 2-4 2-5 2-6 2-7 2-8 2-9 2-10 2-11 2-12	Standard Fluorometer Cuvettes Micro and Semi-Micro Fluorometer Cuvettes Sub-Micro Fluorometer Cuvettes Screw Cap Fluorometer Cuvettes Scattering Tube Cells Fluorometer Flow Cells Fluorometer Micro Flow Cells Water Jacketed Fluorometer Cells Anaerobic Fluorometer Cuvettes Tandem (Divided) Fluorometer Cuvettes Quartz to Glass Grading Fluorometer Cuvettes Cryogenic Fluorometer Cuvettes			
SECTION 3	OTHER CUVETTES			
3-1 3-2 3-3 3-4 3-5	Triangle Cuvettes Liquid Chromatography & Refractometer Flow Cells Dye Laser Cells Laser Absorption Cell Spacers	5		á
SECTION 4	DISPOSABLE/PLASTIC CUVETTES			o,
SECTION 5	ACCESSORIES			Tab
5-1 5-2 5-3 5-4 5-5	Covers, Stoppers, Caps, Special Tops Cell Mounts Cases & Cell Rack Cuvette Washers Mixers/Stirrers			SISTORY
SECTION 6	UV/VIS/NIR Calibration Standards for Spectrophotometers		1	SERVE
6-1 6-2 6-3 6-4 6-5 6-6 6-7	Photometric Accuracy & Absorbance Calibration Wavelength Calibration Stray Light Calibration Spectral Bandwidth Calibration Calibration Kits Fluorescence Reference Standards Microplate Reader Test Plates			S. M.M.M.



About our cuvettes

Each cuvette is thoroughly inspected through our meticulous quality control system, against the industry's highest standards. Rest assured that each cuvette you order has passed this rigorous examination.

Our staff has been supplying reputable researchers and laboratories with standard and specialized cuvettes for many decades.

When you order from FireflySci, Inc., you are ordering from a well established, highly regarded cuvette company. Our cuvettes are 100% guaranteed to meet your quality standards in all regards: Transmission, assembly, matching – everything.

At FireflySci, Inc., quality is number one.

To our customers

We wish to thank you for your support, which clearly indicates your satisfaction with our products. We hope that this new catalog, with its many new, exciting products, will be of help in your professional work.

All our products are guaranteed and are offered at the lowest possible prices without compromising quality.

As one of the many scientists and professionals who now use our products, you know that FireflySci, Inc., specializes in the design and fabrication of the highest quality fire fused cuvettes available on the market today. Our aim is to provide you, our customers, with exactly the type of product your work requires.

FireflySci's scientists, engineers, designers and shop personnel have many years of experience in designing and fabricating cuvettes.

If the exact cell your work requires is not included in this catalog, please contact us. We will also fabricate non-standard cells to your specific requirements.

Your comments and suggestions are always welcome.

Sincerely yours, The Staff at FireflySci, Inc.



Information for manufacturers

The FireflySci factory, has since become well known for its high quality silica glass products, ranging from raw materials to high precision final products. FireflySci has been recognized as a leader in the silica industry for its accumulated abilities, established and innovative technology, and superior quality. We have built reputable quality assurance systems to meet the high standards of our customers, and thusly acquired ISO9001 certification, as well as other domestic quality assurance certifications.

Manufacturers throughout analytical and diagnostic industries, demanding a reliable and technically aggressive silica glass company, have turned to FireflySci for their silica glass components.

Technology innovation and superior quality, combined with excellent service, continues to be the foundation of our philosophy.

Our precision and dedication to quality has allowed us to serve large, world class instrument manufacturers with high caliber optical components for decades.

To see what FireflySci, Inc., can do for your company, please call +1-347-441-4277









Contact Us / Customer Service

We can easily be reached for placing orders, requesting quotations, technical inquiries, customer service, and any other comments or questions you might have.

+1-347-441-4277 +1-347-554-8048 info@fireflysci.com

Mailing Address:

1014 E. 21st Street, Suite A Brooklyn, NY 11210 USA

We can be reached via phone during the following hours: Monday through Friday, 9:00am - 5:00pm
New York time

Ordering Information

Payment:

We accept Visa, MasterCard, and American Express. Please provide number, expiration date, name on card, billing address, and CID number.

We also offer Net 30 Days to established accounts. For an account with us, please send us 4 credit references. We can ship your orders via Net 30 Days once your references have been verified.

Special payment terms apply to international customers. Please inquire.

All prices are FOB/EXW Brooklyn, NY

Shipment, within the USA and Canada:

Our standard mode of shipment is UPS Ground (or UPS Standard to Canada). 1 Day, 2 Day, and 3 Day are also available. Collect shipments to a UPS, TNT, or FedEx account are also available. Collect shipments are shipped insured unless otherwise specified. Terms are FOB Brooklyn, New York.

Shipment, outside the USA and Canada:

Our standard mode of shipment depends on the destination country. We source from a number of carriers to find you the least expensive rates and most reliable service. Standard terms are FCA Brooklyn, New York.

Return policy:

Customers may return cuvettes for refund, exchange, or credit. Returns due to FireflySci error are completely refundable, including shipping.

If items are being returned due to customer error, refund, exchange, or credit can only be issued upon receipt and inspection of cuvettes. The items must be unused to be accepted.

Used cuvettes can also be returned. In these cases, no refund, exchange, or credit is issued until the cuvette has gone through our QC and cleaning departments, at which time a fee is calculated based on the contamination of the cuvette. In some cases, no refund, exchange, nor credit can be provided, and the cuvettes will be shipped back to you as is. Otherwise, a partial credit is applied.

Due to a low volume of returns, no RMA# Is necessary. Please call ahead before you mail the items back.




Material Information

All the cuvettes you order are automatically shipped in optically matched sets. If you want to purchase a new cuvette to match an older one, please let us know.

Also regarding our materials, the following information is available:

- Black quartz information
- Purity
- Chemical compatibility
- Thermal properties
- Fluorescence, refractive index, and other optical qualities

	Usable Ran	ige & Matching Toleran	ices
MATERIAL	ABBREVIATION	USABLE RANGE	MATCHING TOLERANCE
Optical Glass	G	340 - 2,500nm	1% at 350nm
UV Quartz	UV	190 - 2,500nm	1% at 220nm
IR Quartz	IR	220 - 3,500nm	1% at 2,730nm
Pyrex Glass		320 - 2,500nm	
Sapphire		250 - 5,000nm	

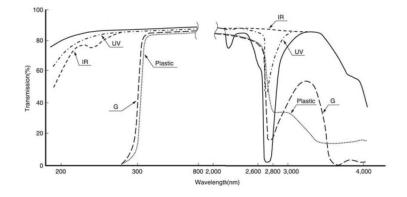

Light Path, Tolerances, and Other Specifications

Regarding general specifications, please note that each part number has its own tolerances. We can also provide information regarding window flatness, polish, and parallelism. Please contact us for this, as well as additional material data listed above.

Frequently Asked Questions

Q: I have a specific instrument, and I need to make sure the cuvettes in this catalog are compatible...

A: The cuvettes in this catalog are guaranteed to work in a myriad of different instruments. Call or email to request a conversion chart between the part numbers of your instruments' cuvettes and our cuvettes.

Q: Is AR coating available?

A. Yes. Please contact us and provide the following information:

- What P/N do you need coated?
 If it is a custom item, please see page 7.
- 2. What wavelength are you working in?
- 3. What is your target percentage of transmission? For example, "98% at 350nm."
- 4. What quantity of the above P/N do you need coated?

Q: Do you manufacture cells for Circular Dichroism (CD)?

A: Yes. Please contact us, and let us know if you need a rectangular or circular cuvettes for your application. We will recommend certain parts to you, and provide more information regarding our cuvettes for CD applications.


Q: What is the Z-dimension (Zd) or Center Height?

A: The Z-dimension is the height from the bottom of the cuvette to the center of its sample chamber. In cases where the cuvette's polished window is very small (where usually the entire cuvette is black aside from the aperture window), the Z-dimension of the cuvette should be the same as the height of the instrument's light beam.

The Z-dimension is usually 8.5 or 15mm, depending on your instrument's manufacturer.

Only certain cuvettes depend on the correct Z-dimension to be used.

These cuvettes have Z-dimension (Zd) information provided in their data table.

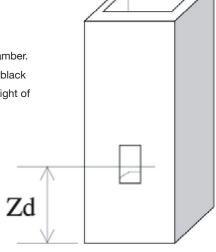

Cleaning and Handling


Cuvettes must be scrupulously clean. Optical surfaces should not be touched as oil smudges are difficult to remove. Purity of solvents used for sample and cleaning is important. Only distilled water or spectrophotometric grade solvents should be used to avoid contamination from impure solvents. As soon as possible after use, cuvettes should be rinsed and soaked in distilled water.

For cleaning, use a mild detergent solution containing no lanolin, oils, or insoluble matter. To remove stubborn contaminants, cuvettes can be soaked for a few minutes in mild sulfuric acid. Never allow cuvettes to soak in hot concentrated acids, alkalis, or other agents which may etch the polished optical surfaces. After cleaning with detergent and/or acid, cuvettes should be thoroughly rinsed with distilled water. Avoid blowing the cuvette with air to dry them. It is better to speed evaporation by means of suction.

Please see page 33 for cuvette cleaners.

Please contact us if you need any further cleaning information.



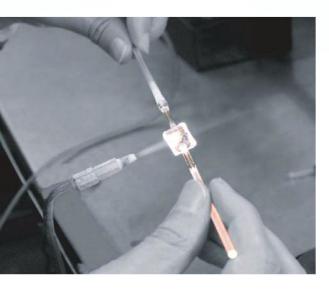

Custom Designs

In order to manufacture your custom quartz glass part, we need the following information:

- 1. Material (Glass, or other (please specify))
- 2. All relevant dimensions and dimension tolerances.


If you do not provide tolerances, our standard tolerances will be used.*

- 3. All critical tolerances and critical areas.*
- 4. Which surfaces to be polished.
- 5. Please fax or email a sketch. A rough drawing is also fine.
- 6. Application information.
- 7. Quantities you need a quotation for.
- 8. Deadline for quotation and/or manufacturing, if applicable.
- *We cannot be responsible for any dimensions not explicitly provided by the customer.







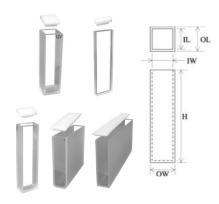

Standard cuvette with PTFE cover

• Two windows polished

Polished = Optical polishing L/P = Light Path

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Materi	al & Part Nur	nbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR
0.5	3.0	12.5	45.5	0.5	10.0	0.220	-	1UV0.5	1IR0.5
1	3.5	12.5	45.5	1.0	10.0	0.400	1G1	1UV1	1IR1
2	4.5	12.5	45.5	2.0	10.0	0.700	1G2	1UV2	1IR2
5	7.5	12.5	45.5	5.0	10.0	1.700	1G5	1UV5	1IR5
10	12.5	12.5	45.5	10.0	10.0	3.500	1G10	1UV10	1IR10
20	22.5	12.5	45.5	20.0	10.0	7.000	1G20	1UV20	1IR20
30	32.5	12.5	45.5	30.0	10.0	10.500	1G30	1UV30	1IR30
40	42.5	12.5	45.5	40.0	10.0	14.000	1G40	1UV40	1IR40
50	52.5	12.5	45.5	50.0	10.0	17.500	1G50	1UV50	1IR50
100	102.5	12.5	45.5	100.0	10.0	35.000	1G100	1UV100	1IR100






TYPE 5

Standard cuvette with PTFE cover, rounded corners

• Two windows polished

	Outer	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	
10	12.5	12.5	45.0	10.0	10.0	3.500	5G10	5UV10	




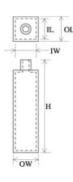

TYPE 11

Standard cuvette with quartz cap

	Outer	Dimensions	(mm)	Inside Dimer	sions (mm)		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
5	7.5	12.5	55.5	5.0	10.0	1.700	11G5	11UV5	11IR5	
10	12.5	12.5	55.5	10.0	10.0	3.500	11G10	11UV10	11IR10	
40	42.5	12.5	55.5	40.0	10.0	14.000	11G40	11UV40	11IR40	
50	52.5	12.5	55.5	50.0	10.0	17.500	11G50	11UV50	11IR50	

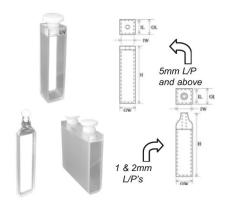


Standard cuvette with PTFE stopper

• Two windows polished

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Materi	al & Part Nui	mbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR
0.5	3.0	12.5	55.5	0.5	10.0	0.220	21G1	21UV1	21IR1
1	3.5	12.5	55.5	1.0	10.0	0.400	21G1	21UV1	21IR1
2	4.5	12.5	55.5	2.0	10.0	0.700	21G2	21UV2	21IR2
5	7.5	12.5	55.5	5.0	10.0	1.700	21G5	21UV5	21IR5
10	12.5	12.5	55.5	10.0	10.0	3.500	21G10	21UV10	21IR10
20	22.5	12.5	55.5	20.0	10.0	7.000	21G20	21UV20	21IR20
30	32.5	12.5	55.5	30.0	10.0	10.500	21G30	21UV30	21IR30
40	42.5	12.5	55.5	40.0	10.0	14.000	21G40	21UV40	21IR40
50	52.5	12.5	55.5	50.0	10.0	17.500	21G50	21UV50	21IR50
100	102.5	12.5	55.5	50.0	10.0	35.500	21G50	21UV50	21IR50



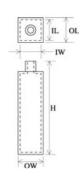


TYPE 31

Standard cuvette with quartz stopper

• Two windows polished

	Outer	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10	12.5	12.5	52.0	10.0	10.0	3.500	31UV10	31IR10	



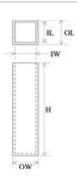


TYPE 508

Tall cuvette with PTFE cover

• Two windows polished

	Outer	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	
10	12.5	12.5	65.0	10.0	10.0	5.500	508G10	508UV10	



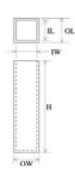


TYPE 509

Tall cuvette with PTFE cover

	Outer	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	
10	12.5	12.5	100	10.0	10.0	10	509G10	509UV10	





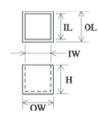

Rectangular absorption cuvette

• Two windows polished

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G
10	13.5	35.0	35.0	10.0	31.5	8.800	523G10

See also Type-37 for a cylindrical cuvette with quartz to glass graded seal tube, page 18.

1-2 Cylindrical Cuvettes



Cylindrical cuvette with PTFE stopper

• Two windows polished




	Oute	r Dimensions (mm)	Inside Dim	ensions (mm)	Material & Part Numbers			
L/P (mm)	LENGTH	DIAMETER	HEIGHT	LENGTH	DIAMETER	CAPACITY (ml)	G	UV	IR
1	3.5	22.0	30.0	1.0	19.0	0.280	32G1	32UV1	32IR1
2	4.5	22.0	30.0	2.0	19.0	0.560	32G2	32UV2	32IR2
5	7.5	22.0	30.0	5.0	19.0	1.400	32G5	32UV5	32IR5
10	12.5	22.0	30.0	10.0	19.0	2.800	32G10	32UV10	32IR10
20	22.5	22.0	30.0	20.0	19.0	5.600	32G20	32UV20	32IR20

TYPE 34

Cylindrical cuvette with 2 PTFE stoppers


• Two windows polished

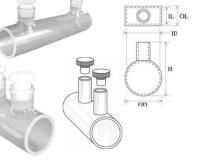

	Oute	r Dimensions (mm)	Inside Dim		Material & Part Numbers			
L/P (mm)	LENGTH	DIAMETER	HEIGHT	LENGTH	DIAMETER	CAPACITY (ml)	G	UV	IR
50	52.5	22.0	35.5	50.0	19.0	14.000	34G50	34UV50	34IR50
100	102.5	22.0	35.5	100.0	19.0	28.000	34G100	34UV100	34IR100
200	202.5	22.0	35.5	200.0	19.0	56.000	34G200	34UV200	34IR200

TYPE 35

Large diameter cylindrical cuvette with PTFE stopper

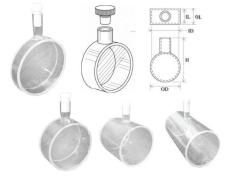

	Oute	Outer Dimensions (mm)			Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	DIAMETER	HEIGHT	LENGTH	DIAMETER	CAPACITY (ml)	G	UV	IR	
10	12.5	50.0	58.0	10.0	47.0	17.000	35G10	35UV10	35IR10	
20	22.5	50.0	58.0	20.0	47.0	34.000	35G20	35UV20	35IR20	
50	52.5	50.0	58.0	50.0	47.0	85.000	35G50	35UV50	35IR50	
100	102.5	50.0	58.0	100.0	47.0	170.000	35G100	35UV100	35IR100	





Short path length cylindrical cuvette with PTFE stopper

• Two windows polished

	Outer Dimensions (mm)		Inside Dim	ensions (mm)	Material & Part Numbers			
L/P (mm)	LENGTH	DIAMETER	HEIGHT	LENGTH	DIAMETER	CAPACITY (ml)	UV	IR
0.05	23.0	22.0	32.0	0.05	13.0	2.150	36UV0.05	36IR0.05
0.1	23.0	22.0	32.0	0.10	13.0	2.150	36UV0.1	36IR0.1
0.2	23.0	22.0	32.0	0.20	13.0	2.180	36UV0.2	36IR0.2
0.5	23.0	22.0	32.0	0.50	13.0	2.220	36UV0.5	36IR0.5
1	23.0	22.0	32.0	1.00	13.0	2.310	36UV1	36IR1
2	23.0	22.0	32.0	2.00	13.0	2.490	36UV2	36IR2
5	23.0	22.0	32.0	5.00	13.0	3.020	36UV5	36IR5




TYPE 521

Long cylindrical cuvette with PTFE stopper

• Two windows polished

	Outer Dimensions (mm)				Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	DIAMETER	HEIGHT	LENGTH	DIAMETER	CAPACITY (ml)	G	UV		
100	102.5	10.0	20.0	100.00	8.0	5.000	521G100	521UV100		

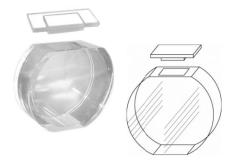


TYPE 522

Cylindrical absorption cuvette with quartz cover

	Outer Dimensions (mm)			Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	DIAMETER	HEIGHT	LENGTH	DIAMETER	CAPACITY (ml)	UV	IR	
10	12.5	35.0	31.5	10.00	31.0	7.000	522UV10	522IR10	





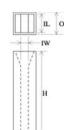

Semi-Micro Cuvette with PTFE cover

• Two windows polished

• Base thickness is 3mm

	Outer Dimensions (mm)		Inside Dimensions (mm)			Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR
5	7.5	12.5	45.0	5.0	4.0	0.700	9G5	9UV5	9IR5
10	12.5	12.5	45.0	10.0	4.0	1.400	9G10	9UV10	9IR10
20	22.5	12.5	45.0	20.0	4.0	2.800	9G20	9UV20	9IR20
30	32.5	12.5	45.0	30.0	4.0	4.200	9G30	9UV30	9IR30
40	42.5	12.5	45.0	40.0	4.0	5.600	9G40	9UV40	9IR40
50	52.5	12.5	45.0	50.0	4.0	7.000	9G50	9UV50	9IR50
100	102.5	12.5	45.0	100.0	4.0	14.000	9G100	9UV100	9IR100



TYPE 9M


Self masking semi-micro cuvette with PTFE cover

• Two windows polished

• Base thickness is 3mm

	Outer Dimensions (mm)			Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR
5	7.5	12.5	45.0	5.0	4.0	0.700	9MG5	9MUV5	9MIR5
10	12.5	12.5	45.0	10.0	4.0	1.400	9MG10	9MUV10	9MIR10
20	22.5	12.5	45.0	20.0	4.0	2.800	9MG20	9MUV20	9MIR20
40	42.5	12.5	45.0	40.0	4.0	5.600	9MG40	9MUV40	9MIR40
50	52.5	12.5	45.0	50.0	4.0	7.000	9MG50	9MUV50	9MIR50
100	102.5	12.5	45.0	100.0	4.0	14.000	9MG100	9MUV100	9MIR100



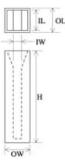


TYPE 9B

Variable capacity micro/semi-micro cuvette with PTFE cover

- Two windows polished
- Capacity changes based on inside width
- Base thickness is 9mm
- No black masking

	Outer	Outer Dimensions (mm)			Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR		
10	12.5	12.5	45.0	10.0	1.0	0.4	9BUV10X1	9BIR10X1		
10	12.5	12.5	45.0	10.0	3.0	1.1	9BUV10X3	9BIR10X3		
10	12.5	12.5	45.0	10.0	5.0	1.8	9BUV10X5	9BIR10X5		



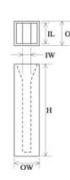

TYPE 9BM

Variable capacity self masking micro/semi-micro cuvette with PTFE cover


• Two windows polished • Capacity changes based on inside width • Base thickness is 9mm

	Outer	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10	12.5	12.5	45.0	10.0	1.0	0.4	9BMUV10X1	9BMIR10X1	
10	12.5	12.5	45.0	10.0	3.0	1.1	9BMUV10X3	9BMIR10X3	
10	12.5	12.5	45.0	10.0	5.0	1.8	9BMUV10X5	9BMIR10X5	





1-3

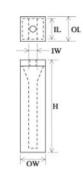

TYPE 29

Semi-micro cuvette with PTFE stopper

- Two windows polished
- Base thickness is 3mm

	Outer Dimensions (mm)			Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR
5	7.5	12.5	49.0	5.0	4.0	0.700	29G5	29UV5	29IR5
10	12.5	12.5	49.0	10.0	4.0	1.400	29G10	29UV10	29IR10
20	22.5	12.5	49.0	20.0	4.0	2.800	29G20	29UV20	29IR20
30	32.5	12.5	49.0	30.0	4.0	4.200	29G30	29UV30	29IR30
40	42.5	12.5	49.0	40.0	4.0	5.600	29G40	29UV40	29IR40
50	52.5	12.5	49.0	50.0	4.0	7.000	29G50	29UV50	29IR50




TYPE 29M

IL OL

Self masking semi-micro cuvette with PTFE stopper

- Two windows polished
- Base thickness is 3mm

	Outer Dimensions (mm)			Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR
5	7.5	12.5	49.0	5.0	4.0	0.700	29MG5	29MUV5	29MIR5
10	12.5	12.5	49.0	10.0	4.0	1.400	29MG10	29MUV10	29MIR10
20	22.5	12.5	49.0	20.0	4.0	2.800	29MG20	29MUV20	29MIR20
40	42.5	12.5	49.0	40.0	4.0	5.600	29MG40	29MUV40	29MIR40
50	52.5	12.5	49.0	50.0	4.0	7.000	29MG50	29MUV50	29MIR50

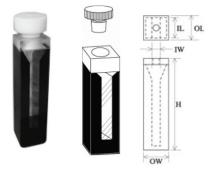

TYPE 29BM

Semi-Micro Cuvette with PTFE cover

- Two windows polished
- Base thickness is 9mm

	Outer	Outer Dimensions (mm)			sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
5	7.5	12.5	49.0	5.0	4.0	0.600	29BMUV5	29BMIR5	
10	12.5	12.5	49.0	10.0	4.0	1.200	29BMUV10	29BMIR10	
20	22.5	12.5	49.0	20.0	4.0	2.400	29BMUV20	29BMIR20	
30	32.5	12.5	49.0	30.0	4.0	3.600	29BMUV30	29BMIR30	
40	42.5	12.5	49.0	40.0	4.0	4.800	29BMUV40	29BMIR40	



Multi-Micro Cell Array

- 8-chambered multi-micro array of cells for doing multiple measurements.
- Typically used in Shimadzu machines.

	Outer Dimensions (mm)			Inside Dimensions (mm)			Material & Part Numbers	
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	
10mm	12.5mm	74.5mm	9.0mm	10.0	10.0	0.120	8UV10	




TYPE 16

Multi-Micro Cell Array

• 16-chambered multi-micro array of cells for doing multiple measurements.

	Outer	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers	
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	
5	12.5	74.5	9.0	5.0	2.0	0.50	16UV5	
10	12.5	74.5	9.0	10.0	2.0	0.100	16UV10	



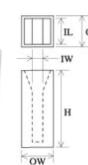

TYPE 17

Short micro cuvette with PTFE cover

• Two windows polished


Base thickness is 3mm

• Height is 25mm

	Outer	Outer Dimensions (mm)			sions (mm)	Material & Part Numbers				
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
5	7.5	12.5	25.0	5.0	2.0	0.200	17G5	17UV5	17IR5	
10	12.5	12.5	25.0	10.0	2.0	0.400	17G10	17UV10	17IR10	
20	22.5	12.5	25.0	20.0	2.0	0.800	17G20	17UV20	17IR20	
30	32.5	12.5	25.0	30.0	2.0	1.200	17G30	17UV30	17IR30	
40	42.5	12.5	25.0	40.0	2.0	1.600	17G40	17UV40	17IR40	
50	52.5	12.5	25.0	50.0	2.0	2.000	17G50	17UV50	17IR50	





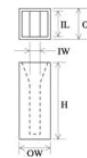

TYPE 17M

or. Self masking short micro cuvette with PTFE cover


• Two windows polished

• Base thickness is 3mm

• Height is 25mm

	Outer Dimensions (mm)			Inside Dimensions (mm)			Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
5	7.5	12.5	25.0	5.0	2.0	0.200	17MG5	17MUV5	17MIR5	
10	12.5	12.5	25.0	10.0	2.0	0.400	17MG10	17MUV10	17MIR10	
20	22.5	12.5	25.0	20.0	2.0	0.800	17MG20	17MUV20	17MIR20	
30	32.5	12.5	25.0	30.0	2.0	1.200	17MG30	17MUV30	17MIR30	
40	42.5	12.5	25.0	40.0	2.0	1.600	17MG40	17MUV40	17MIR40	
50	52.5	12.5	25.0	50.0	2.0	2.000	17MG50	17MUV50	17MIR50	

IL OL

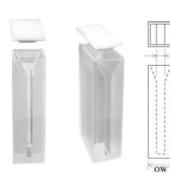

1-3

TYPE 18


Micro cuvette with PTFE cover

• Two windows polished

• Base thickness is 3mm

	Outer Dimensions (mm)			Inside Dimer		Material & Part Numbers				
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
5	7.5	12.5	45.0	5.0	2.0	0.350	18G5	18UV5	18IR5	
10	12.5	12.5	45.0	10.0	2.0	0.700	18G10	18UV10	18IR10	
20	22.5	12.5	45.0	20.0	2.0	1.400	18G20	18UV20	18IR20	
30	32.5	12.5	45.0	30.0	2.0	2.100	18G30	18UV30	18IR30	
40	42.5	12.5	45.0	40.0	2.0	2.800	18G40	18UV40	18IR40	
50	52.5	12.5	45.0	50.0	2.0	3.500	18G50	18UV50	18IR50	
100	102.5	12.5	45.0	100.0	2.0	7.000	18G100	18UV100	18IR100	



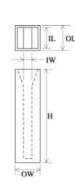

TYPE 18M


Self masking micro-cuvette with PTFE cover

• Two windows polished

• Base thickness is 3mm

	Outer I	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
5	7.5	12.5	45.0	5.0	2.0	0.350	18MG5	18MUV5	18MIR5	
10	12.5	12.5	45.0	10.0	2.0	0.700	18MG10	18MUV10	18MIR10	
20	22.5	12.5	45.0	20.0	2.0	1.400	18MG20	18MUV20	18MIR20	
40	42.5	12.5	45.0	40.0	2.0	2.800	18MG40	18MUV40	18MIR40	
50	52.5	12.5	45.0	50.0	2.0	3.500	18MG50	18MUV50	18MIR50	



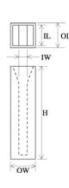

TYPE 18B


Bottom raised micro cuvette with PTFE cover


• Two windows polished

• Base thickness is 9mm

No black masking

	Outer	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
5	7.5	12.5	45.0	5.0	2.0	0.280	18BUV5	18BIR5	
10	12.5	12.5	45.0	10.0	2.0	0.560	18BUV10	18BIR10	
20	22.5	12.5	45.0	20.0	2.0	1.120	18BUV20	18BIR20	
40	42.5	12.5	45.0	40.0	2.0	2.240	18BUV40	18BIR40	
50	52.5	12.5	45.0	50.0	2.0	2.800	18BUV50	18BIR50	

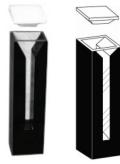

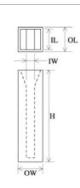

TYPE 18BM


Bottom raised self masking micro cuvette with PTFE cover

• Two windows polished

Base thickness is 9mm

	Outer	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
5	7.5	12.5	45.0	5.0	2.0	0.280	18BMUV5	18BMIR5	
10	12.5	12.5	45.0	10.0	2.0	0.560	18BMUV10	18BMIR10	
20	22.5	12.5	45.0	20.0	2.0	1.120	18BMUV20	18BMIR20	



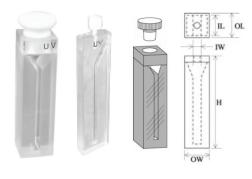


Micro cuvette with PTFE stopper

• Two windows polished

• Base thickness is 3mm

	Outer Dimensions (mm)			Inside Dimensions (mm)			Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
1	3.5	12.5	49.0	1.0	2.0	0.070	30G1	30UV1	30IR1	
2	4.5	12.5	49.0	2.0	2.0	0.140	30G2	30UV2	30IR2	
5	7.5	12.5	49.0	5.0	2.0	0.350	30G5	30UV5	30IR5	
10	12.5	12.5	49.0	10.0	2.0	0.700	30G10	30UV10	30IR10	
30	32.5	12.5	49.0	30.0	2.0	2.100	30G30	30UV30	30IR30	
40	42.5	12.5	49.0	40.0	2.0	2.800	30G40	30UV40	30IR40	
50	52.5	12.5	49.0	50.0	2.0	3.500	30G50	30UV50	30IR50	

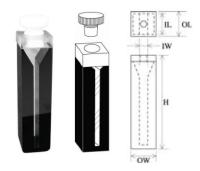

TYPE 30M


Self masking micro cuvette with PTFE stopper

• Two windows polished

• Base thickness is 3mm

	Outer	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
5	7.5	12.5	49.0	5.0	2.0	0.350	30MG5	30MUV5	30MIR5	
10	12.5	12.5	49.0	10.0	2.0	0.700	30MG10	30MUV10	30MIR10	
20	22.5	12.5	49.0	20.0	2.0	1.400	30MG20	30MUV20	30MIR20	
40	42.5	12.5	49.0	40.0	2.0	2.800	30MG40	30MUV40	30MIR40	



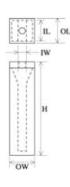

TYPE 30BM

Bottom raised self masking micro cuvette with PTFE stopper

• Two windows polished

Base thickness is 9mm

	Outer	Dimensions	(mm)	Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
5	7.5	12.5	49.0	5.0	2.0	0.300	30BMUV5	30BMIR5	
10	12.5	12.5	49.0	10.0	2.0	0.600	30BMUV10	30BMIR10	
20	22.5	12.5	49.0	20.0	2.0	1.200	30BMUV20	30BMIR20	



TYPE 520M

Short self masking micro cuvette with PTFE stopper

• Two windows polished

• Base thickness is 3mm

	Outer Dimensions (mm)			Inside Dimen	sions (mm)	Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
5	7.5	12.5	29.0	5.0	2.0	0.200	520MUV5	520MIR5	
10	12.5	12.5	29.0	10.0	2.0	0.400	520MUV10	520MIR10	
20	22.5	12.5	29.0	20.0	2.0	0.800	520MUV20	520MIR20	





TYPE 701M

Sub-Micro Cuvettes

Sub-micro self masking cuvette with PTFE threaded stopper

• Two windows polished

	Outer I	Dimensions	(mm)	Chambe	r Dimensio	ns (mm)			Material & Part Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (μl)	UV
10	12.5	12.5	45.0	10	1.0	1.0	8.5	10	701MUV10.10A
10	12.5	12.5	45.0	10	1.0	1.0	15.0	10	701MUV10.10B
10	12.5	12.5	45.0	10	2.0	2.5	20.0	10	701MUV10.10C
10	12.5	12.5	45.0	10	2.0	2.5	8.5	50	701MUV10.50A
10	12.5	12.5	45.0	10	2.0	5.0	15.0	50	701MUV10.50B
10	12.5	12.5	45.0	10	2.0	5.0	20.0	50	701MUV10.50C
10	12.5	12.5	45.0	10	2.0	8.0	8.5	100	701MUV10.100A
10	12.5	12.5	45.0	10	2.0	8.0	15.0	100	701MUV10.100B
10	12.5	12.5	45.0	10	2.0	8.0	20.0	100	701MUV10.100C
10	12.5	12.5	45.0	10	2.0	8.0	8.5	160	701MUV10.160A
10	12.5	12.5	45.0	10	2.0	8.0	15.0	160	701MUV10.160B
10	12.5	12.5	45.0	10	2.0	8.0	20.0	160	701MUV10.160C







1-5 Ultra Micro Cuvette

Micro Focus Cell

Ultra-micro low volume cell with a 1.9µl capacity

	Outer	Dimensions	s (mm)	Chamber Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH WIDTH HEIGHT		LENGTH	ZD (mm)	CAPACITY (μΙ)	UV			
1	12.5	12.5	45.0	1.0	8.5- 20.0	1.9	MFC-1		

INCLUDED IN THE KIT IS THE FOLLOWING:

- A 12.5x12.5mm micro focus adapter that fits into any spectrophotometer.
- Cleaning swabs.

• A UV Quartz micro tray to inject a sample in.

- Instructions.
- Spacers to adjust the micro focus adapter to the proper Z-Dimension in your spectrophotometer.

HERE IS A BASIC GUIDE TO USING THE MICRO FOCUS CELL CONVERSION KIT:

1. To use the conversion kit all you have to do is check the Z-dimension of your machine.

The Z-Dimension of the Micro Focus adapter is 8.5 mm so if you have a machine with an 8.5 or 20mm z-dimension you'll need to use the appropriate spacer.

2. Next insert the Micro Focus Adapter into the spectrophotometer. Pull out the UV Quartz micro chamber.

The micro chamber is a 1 mm lightpath cuvette that holds a sample of 1.9µl.

- 3. Inject your sample into the quartz chamber and insert this into the Micro Focus Adapter.
- 4. Run your experiment.
- 5. Celebrate how easy that was to do and how much money you saved.

APPLICATIONS FOR THE MICRO FOCUS CELL:

- DNA Measurement
- RNA Measurement
- DNA Quantification
- Nucleic Acid Analysis
- Protein Concentration
- All small volume applications (0.3-1.9µl)
- Only the Micro Focus Cell can offer these guarantees:
- Quick: Load and empty your sample into the cell without even removing it from the machine.

This equals consistent alignment through multiple tests.

• Guaranteed to work in your spectrophotometer (any Z dimension/center height is compatible) Leak proof design. Our job is design and manufacture the most innovative and quality cuvettes. Your job is to use them and change the world: 0.1µl at a time. Let us know how else we can help.





TYPE 34S

Large Cylindrical Polarimeter Cell with Screw Caps

• Includes 2 SC2 open screw caps and two SC1 closed screw caps

	Oute	r Dimension	ns (mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	HEIGHT	DIAMETER	DIAMETER	LENGTH	CAPACITY (ml)	UV		
50	52.5	35.5	22.0	19.0	50.0	14.000	34SUV50		
100	102.5	35.5	22.0	19.0	100.0	28.000	34SUV100		





TYPE 41

Standard cuvette with screw cap


- Two windows polished
- Height includes threading, but not screw cap
- Cuvette comes with both open and closed screw caps
- Each cap contains a PTFE septa seal
- Suitable for anaerobic work

	Outer Dimensions (mm)			Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)			LENGTH	WIDTH	CAPACITY (ml)	G	UV		
1	3.5	12.5	55.0	1.0	10.0	0.400	-	41UV1	
2	4.5	12.5	55.0	2.0	10.0	0.700	-	41UV2	
5	7.5	12.5	55.0	5.0	10.0	1.700	41G5	41UV5	
10	12.5	12.5	55.0	10.0	10.0	3.500	41G10	41UV10	
50	52.5	12.5	55.0	50.0	10.0	17.500	41G50	41UV50	
100	00 102.5 12.5 55.0	100.0	10.0	35.000	41G100	41UV100			







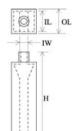

TYPE 46

Variable capacity micro/semi-micro cuvette with screw cap

- Two windows polished
- Height includes threading, but not screw cap
- Cuvette comes with both open and closed screw caps
- Each cap contains a PTFE septa seal
- Suitable for anaerobic work
- Base thickness is 3mm

	Outer	Dimensions	(mm)	Inside Dimer	nsions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH WIDTH		CAPACITY (ml)	UV		
10	12.5	12.5	55.0	10.0	1.0	0.350	46UV10X1		
10	12.5	12.5	55.0	10.0	2.0	0.700	46UV10X2		
10	12.5	12.5	55.0	10.0	3.0	1.050	46UV10X3		
10	12.5	12.5	55.0	10.0	4.0	1.400	46UV10X4		
10	12.5	12.5	55.0	10.0	5.0	1.750	46UV10X5		

Flow Cells/TD (Tablet Dissolution) Cells



TYPE 44

Flow cell with detachable inlet/outlet tubes, 1.6ml

- Two windows polished
- Inlet/outlet tubes: OD4 x ID2mm
- Height is without tubes
- Tubes are made of quartz

	Outer D	Dimensions	s (mm)	Chambe	r Dimensio	ns (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR
10	12.5	12.5	50.0	10.0	4.0	40.0	1.600	44UV10	44IR10



TYPE 45

Flow cell with detachable inlet/outlet tubes, 2.8ml

- Two windows polished
- Inlet/outlet tubes: OD4 x ID2mm
- Height is without tubes
- Tubes are made of quartz

	Outer Dimensions (mm)			Chambe	r Dimensio	ons (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR	
10	12.5	12.5	50.0	10.0	7.0	40.0	2.800	45UV10	45IR10	



TYPE 48

Flow-through cuvette with side tubes

- Two windows polished
- Inlet/outlet tubes: OD5 x ID3 x L15mm
- Length is without tubes
- Tubes are made of quartz

	Outer E	Dimension	s (mm)	Chambe	r Dimensio	ons (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR	
0.05	2.6	12.5	45.0	0.05	9.0	38.0	0.015	48UV0.05	48IR0.05	
0.1	2.6	12.5	45.0	0.1	9.0	38.0	0.030	48UV0.1	48IR0.1	
0.5	3.0	12.5	45.0	0.5	9.0	38.0	0.150	48UV0.5	48IR0.5	
1	3.5	12.5	45.0	1.0	9.0	38.0	0.300	48UV1	48IR1	
2	4.5	12.5	45.0	2.0	9.0	38.0	0.600	48UV2	48IR2	
5	7.5	12.5	45.0	5.0	9.0	38.0	1.500	48UV5	48IR5	
10	12.5	12.5	45.0	10.0	9.0	38.0	3.000	48UV10	48IR10	



TYPE 58

Flow-through cell with top tubes

- Two windows polished
- Inlet/outlet tubes: OD4 x ID2 x L10mm
- Height is without tubes
- Tubes are made of quartz

	Outer Dimensions (mm)			Chambe	r Dimensio	ons (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR	
10	12.5	12.5	45.0	10.0	7.0	36.0	2.520	58UV10	58IR10	





Flow-through cell with top and bottom tubes

- Two windows polished
- Inlet/outlet tubes: OD5 x ID3 x L10mm
- Height is without tubes
- Tubes are made of quartz

	Outer E	Dimension	s (mm)	Chamber	Dimensio	ns (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR
5	12.5	12.5	65.0	5.0	10.0	35.0	1.850	501UV5	501IR5
10	12.5	12.5	45.0	10.0	10.0	33.0	3.300	501UV10	501UV10

1-8 Micro & Self Masking Flow Cells



Micro flow through cell, 2mm L/P

- Two windows polished
- PTFE inlet/outlet tubes: OD2 x ID1mm
- Height is without tubes

	Outer D	Dimensions	s (mm)	Chambe	r Dimensio	ns (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR	
2	3.0	3.0	7.6	2.0	2.0	7.6	0.030	74UV2	74IR2	



TYPE 513M

Micro flow through cell, 2mm L/P

- Two windows polished PTFE inlet/outlet tubes: OD4 x ID2mm
- Height is without tubes

	Outer l	Dimensions	(mm)	Chamber Dime	ensions (mm)	Material & Part Numbers			
L/P (mm)	LENGTH WIDTH HEIGHT		LENGTH WIDTH		ZD (mm)	CAPACITY (ml)	UV	IR	
10	12.5	12.5	38.0	10.0	3.0	8.5	0.070	513MUV10	513MIR10



TYPE 59M

Self masking micro flow cell

- Two windows polished
- Inlet/outlet tubes: OD3 x ID1.5 x L8mm
- Height is without tubes
- Tubes are made of quartz

		Outer I	Dimension	s (mm)	Chambei	r Dimensio	ons (mm)	Material & Part Number			
	L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (ml)	UV	IR
Ī	10	12.5	12.5	45.0	10.0	4.0	11.0	15.0	0.440	59MUV10	59MIR10





TYPE 500M

Fully masked ultra-micro flow cell with cylindrical chamber

- Two windows polished
- Inlet/outlet tubes: OD2 x ID1 x L5mm
- Height is without inlet/outlet tubes

	Outer	Dimensions	(mm)	Chamber Dime	ensions (mm)				laterial & t Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	ZD (mm)	CAPACITY (ml)	UV	IR
10	12.5	12.5	38.0	10.0	3.0	8.5	0.070	500MUV10	500MIR10

1-9 **Water Jacketed Cuvettes**



TYPE 50

Water jacketed cylindrical cuvette with PTFE stopper

- Two windows polished
- Tubes to water jacket: OD5 x ID3 x L25mm
- Height is without inlet/outlet tubes

	Outer	Dimensions	(mm)	Inside Dimer	sions (mm)	Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
5	7.5	23.0	32.0	5.0	13.0	0.660	50UV5	50IR5	
10	12.5	23.0	32.0	10.0	13.0	1.320	50UV10	50IR10	
50	52.5	23.0	32.0	50.0	13.0	6.600	50UV50	50IR50	



TYPE 51H

Water jacketed cuvette with horizontal tubes

- Two windows polished
- Tubes to water jacket: OD4 x ID2 x L12.5mm
- Width is without tubes

	Outer E	Dimensions	s (mm)	Chambe	r Dimensio	ons (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR
10	12.5	12.5	50.0	10.0	4.0	45.0	1.500	51HUV10	51HIR10



TYPE 51V

Water jacketed cuvette with vertical tubes

- Two windows polished
- Tubes to water jacket: OD4 x ID2 x L12.5mm • Height is without tubes

	Outer E	Dimensions	s (mm)	Chambei	Dimensio	ns (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR
10	12.5	12.5	50.0	10.0	4.0	45.0	1.500	51VUV10	51VIR10





Water jacketed cuvette with PTFE stopper, 2.4ml

- Two windows polished
- Tubes to water jacket: OD4 x ID2 x L12.5mm
- Width is without tubes
- Height is without stopper

	Outer D	Dimensions	s (mm)	Chambe	r Dimensio	ns (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR
10	12.5	12.5	50.0	10.0	4.0	42.0	1.400	53UV10	53IR10





TYPE 54

Water jacketed cuvette with PTFE stopper, 1.2ml

- Two windows polished
- Tubes to water jacket: OD4 x ID2 x L10mm
- Length is without tubes
- Height is without stopper

	Outer D	Dimension	s (mm)	Chambe	r Dimensio	ns (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR
10	12.5	12.5	48.0	10.0	4.0	36.0	1.200	54UV10	54IR10

1-10 Anaerobic Cuvettes



TYPE 26

Short anaerobic cuvette with glass pouch

- Two windows polished
- Overall height is approximately 65mm
- Height is without glass pouch

	Outer E	Dimensions	s (mm)	Chamber	Dimensio	ns (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR	
2	4.5	12.5	45.0	2.0	10.0	30.5	0.500	26UV2	26IR2	
5	7.5	12.5	45.0	5.0	10.0	30.5	1.250	26UV5	26IR5	
10	12.5	12.5	45.0	10.0	10.0	30.5	2.500	26UV10	26IR10	



TYPE 28

Standard anaerobic cuvette with glass pouch

- Two windows polished
- Overall height is approximately 125mm
- Height is without glass pouch

	Outer E	Dimensions	s (mm)	Chambei	Dimensio	ons (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR	
2	4.5	12.5	84.0	2.0	10.0	44.0	0.700	28UV2	28IR2	
5	7.5	12.5	84.0	5.0	10.0	44.0	1.750	28UV5	28IR5	
10	12.5	12.5	84.0	10.0	10.0	44.0	3.500	28UV10	28IR10	





Anaerobic Cuvettes

Bubble forming anaerobic cuvette with glass pouch

- Two windows polished
- Overall height is approximately 102mm
- Height is without glass pouch

	Outer E	Dimensions	s (mm)	Chambe	r Dimensio	ns (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR
10	12.5	12.5	76.5	10	10	45.0	3.500	71UV10	71IR10

1-11 **Tandem (Divided) Cuvettes**



TYPE 25A

Tandem cuvette

- Two windows & partition polished
- Height of partition is approximately 37mm

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)	Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10	12.5	12.5	45.0	2x4.5	10	2 x 1.300	25AUV10	25AIR10	

1-12 **Demountable Cuvettes**

*A20 cell mount not included. See page 48

TYPE 19

U-demountable cuvette with removable window

	Outer E	Dimensions	s (mm)	Chambe	Dimensio	ons (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR
0.05	2.55	12.5	45.0	0.05	10.0	41.5	0.020	19UV0.05	19IR0.05
0.1	2.6	12.5	45.0	0.1	10.0	41.5	0.040	19UV0.1	19IR0.1
0.2	2.7	12.5	45.0	0.2	10.0	41.5	0.080	19UV0.2	19IR0.2
0.5	3.0	12.5	45.0	0.5	10.0	41.5	0.200	19UV0.5	19IR0.5
1	3.5	12.5	45.0	1.0	10.0	41.5	0.400	19UV1	19IR1
2	4.5	12.5	45.0	2.0	10.0	41.5	0.720	19UV2	19IR2





*A20 cell mount not included. See page 48

O-demountable cuvette with removable window

• Two windows polished


	Outer E	Dimensions	s (mm)	Chamber	r Dimensio	ns (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR	
0.05	2.55	12.5	45.0	0.05	10.0	38.0	0.018	20UV0.05	20IR0.05	
0.1	2.6	12.5	45.0	0.1	10.0	38.0	0.036	20UV0.1	20IR0.1	
0.2	2.7	12.5	45.0	0.2	10.0	38.0	0.072	20UV0.2	20IR0.2	
0.5	3.0	12.5	45.0	0.5	10.0	38.0	0.180	20UV0.5	20IR0.5	
1	3.5	12.5	45.0	1.0	10.0	38.0	0.360	20UV1	20IR1	
2	4.5	12.5	45.0	2.0	10.0	38.0	0.720	20UV2	20IR2	



- Two windows optically polished
- Disc-shaped demountable spectrophotometer cuvette with removable window

	Outer Dimen	sions (mm)	Chamber Dimensions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH WIDTH		WIDTH	CAPACITY (ml)	UV		
0.1	2.6	22.0	16	0.02	42UV0.1		
0.2	2.7	22.0	16	0.04	42UV0.2		
0.5	3.0	22.0	16	0.08	42UV0.5		
1	3.5	22.0	16	0.16	42UV1		
2	4.5	22.0	16	0.32	42UV2		



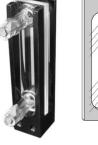

*A20 cell mount not included. See page 48

TYPE 49

O-demountable cuvette with removable window

- Two windows polished
- Inlet/outlet tubes: OD5 x ID3 x L15mm
- Tubes are made of quartz Outside length does not include tubes

	Outer E	Dimensions	s (mm)	Chambei	r Dimensio	ons (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR
0.1	2.6	12.5	45.0	0.1	9.0	38.0	0.030	49UV0.1	49IR0.1
0.2	2.7	12.5	45.0	0.2	9.0	38.0	0.060	49UV0.2	49IR0.2
0.5	3.0	12.5	45.0	0.5	9.0	38.0	0.150	49UV0.5	49IR0.5
1	3.5	12.5	45.0	1.0	9.0	38.0	0.300	49UV1	49IR1
2	4.5	12.5	45.0	2.0	9.0	38.0	0.600	49UV2	49IR2



1-13

TYPE 37

Micro cylindrical cuvette with graded seal tube

- Two windows polished
- Quartz or glass to glass graded seal tube: OD5 x ID3mm
- Overall height is approximately 92mm

	Outer Dime	nsions (mm)	Inside Dim	ensions (mm)	Material & Part Numbers			
L/P (mm)	LENGTH DIAMETER		LENGTH	DIAMETER	CAPACITY (ml)	G	UV	IR
1	3.5	22.0	1.0	19.0	0.280	37G1	37UV1	37IR1
2	4.5	22.0	2.0	19.0	0.570	37G2	37UV2	37IR2
5	7.5	22.0	5.0	19.0	1.420	37G5	37UV5	37IR5
10	12.5	22.0	10.0	19.0	2.830	37G10	37UV10	37IR10
20	22.5	22.0	20.0	19.0	5.670	37G20	37UV20	37IR20
50	52.5	22.0	50.0	19.0	14.000	37G50	37UV50	37IR50
100	102.5	22.0	100.0	19.0	28.000	37G100	37UV100	37IR100



TYPE 61

Standard cuvette with graded seal tube

- Two windows polished
- Quartz to glass graded seal tube: OD8 x ID6mm
- Overall height is approximately 125mm

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR		
0.5	3.0	12.5	45.0	0.5	10.0	0.200	61UV0.5	61IR0.5		
1	3.5	12.5	45.0	1.0	10.0	0.400	61UV1	61IR1		
2	4.5	12.5	45.0	2.0	10.0	0.700	61UV2	61IR2		
5	7.5	12.5	45.0	5.0	10.0	1.700	61UV5	61IR5		
10	12.5	12.5	45.0	10.0	10.0	3.500	61UV10	61IR10		



TYPE 504

Gas absorption cuvette with graded seal tube

- Two windows polished
- Quartz to glass graded seal tube: OD8 x ID6mm
- Overall height is approximately 84mm

	Outer Dime	ensions (mm)	Inside Dim	ensions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	DIAMETER	LENGTH	DIAMETER	CAPACITY (ml)	UV	IR	
74	80.0	22.0	74.0	19.0	21.000	504UV74	504IR74	



TYPE 505

Semi-micro cuvette with graded seal tube

- Two windows polished
- Quartz to glass graded seal tube: OD8 x ID6mm
- Base thickness is 3mm
- Overall height is approximately 125mm

	Outer	Dimensions	(mm)	Inside Dimer	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10	12.5	12.5	45.0	10.0	4.0	1.400	505UV10	505IR10	



Spectrophotometer Cuvettes



TYPE 64

Cryogenic cuvette with PTFE cover

- Two windows polished
- All joints are excessively fused to round corners
- Minimum working area of windows is 4mm wide

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10	12.5	12.5	45.0	10.0	10.0	3.500	64UV10	64IR10	



TYPE 65

Cryogenic cuvette with PTFE stopper

- Two windows polished
- All joints are excessively fused to round corners
- Minimum working area of windows is 4mm wide

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
2	4.5	12.5	45.0	2.0	10.0	0.700	65UV2	65IR2	
10	12.5	12.5	45.0	10.0	10.0	3.500	65UV10	65IR10	

1-15 Tablet Dissolution Cells



TYPE 601

Tablet Dissolution Cells 10 - 100µl

- Two windows polished
- With M6 threaded connectors and PTFE tubing
- In metal housing
- All black quartz; Size of polished area is 2 x 5mm

	Outer I	Dimension	s (mm)	Chamber	r Dimensio	ns (mm)			Part Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (ml)	UV
1	12.5	12.5	35.0	1.0	2.0	5.0	8.5	10	601UV1A
1	12.5	12.5	35.0	1.0	2.0	5.0	15.0	10	601UV1B
2	12.5	12.5	35.0	2.0	2.0	5.0	8.5	20	601UV2A
2	12.5	12.5	35.0	2.0	2.0	5.0	15.0	20	601UV2B
5	12.5	12.5	35.0	5.0	2.0	5.0	8.5	50	601UV5A
5	12.5	12.5	35.0	5.0	2.0	5.0	15.0	50	601UV5B
10	12.5	12.5	35.0	10.0	2.0	5.0	8.5	100	601UV10A
10	12.5	12.5	35.0	10.0	2.0	5.0	15.0	100	601UV10B



TYPE 602

Tablet Dissolution Cells 160µl

- Two windows polished
- With M6 threaded connectors and PTFE tubing
- In metal housing
- All black quartz; Size of polished area is 2 x 8mm

	Outer I	Dimension	s (mm)	Chamber	Dimensio	Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (ml)	UV
10	12.5	12.5	35.0	10.0	2.0	8.0	8.5	160	602UV10A
10	12.5	12.5	35.0	10.0	2.0	8.0	15.0	160	602UV10B



1-15



Tablet Dissolution Cells

TYPE 604

Tablet Dissolution Cells 5 ~ 100mm Light Paths

- Two windows polished
- All black quartz; Size of polished area is 3.5 x 11mm
- Tubes are made of quartz, no threading

	Outer I	Dimension	s (mm)	Chamber	Dimensio	ns (mm)		Material & Part Numbers	
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (µI)	UV
50	52.5	12.5	35.0	50.0	3.5	11.0	8.5	195	604UV50A
100	102.5	12.5	35.0	100.0	3.5	11.0	8.5	390	604UV100A



TYPE 605

Tablet Dissolution Cells 62 & 124µl

- Two windows polished
- With M6 threaded connectors and PTFE tubing
- In metal housing
- All black quartz; Size of polished area is 3.5 x 17.5mm

	Outer I	Dimension	s (mm)	Chamber	r Dimensio	ns (mm)			Material & Part Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (μl)	UV
0.1	12.5	12.5	35.0	0.1	3.5	17.5	8.5	5	605UV0.1A
0.1	12.5	12.5	35.0	0.1	3.5	17.5	15	5	605UV0.1B
0.1	12.5	12.5	35.0	0.1	3.5	17.5	20	5	605UV0.1C
0.2	12.5	12.5	35.0	0.2	3.5	17.5	8.5	12	605UV0.2A
0.2	12.5	12.5	35.0	0.2	3.5	17.5	15	12	605UV0.2B
0.2	12.5	12.5	35.0	0.2	3.5	17.5	20	12	605UV0.2C
0.5	12.5	12.5	35.0	0.5	3.5	17.5	8.5	25	605UV0.5A
0.5	12.5	12.5	35.0	0.5	3.5	17.5	15	25	605UV0.5B
0.5	12.5	12.5	35.0	0.5	3.5	17.5	20	25	605UV0.5C
1	12.5	12.5	35.0	1.0	3.5	17.5	8.5	62	605UV1A
1	12.5	12.5	35.0	1.0	3.5	17.5	15	62	605UV1B
1	12.5	12.5	35.0	1.0	3.5	17.5	20	62	605UV1C
2	12.5	12.5	35.0	2.0	3.5	17.5	8.5	124	605UV2A
2	12.5	12.5	35.0	2.0	3.5	17.5	15	124	605UV2B
2	12.5	12.5	35.0	2.0	3.5	17.5	20	124	605UV2C



TYPE 606

Tablet Dissolution Cells 195 & 390µl

- Two windows polished
- With M6 threaded connectors and PTFE tubing
- In metal housing
- All black quartz; Size of polished area is 3.5 x 11mm

	Outer Dimensions (mm)			Chamber	Dimensio	ns (mm)			Part Numbers	
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (μl)	UV	
5	12.5	12.5	35.0	5.0	3.5	11.0	8.5	195	606UV5A	
5	12.5	12.5	35.0	5.0	3.5	11.0	15.0	195	606UV5B	
10	12.5	12.5	35.0	10.0	3.5	11.0	8.5	390	606UV10A	
10	12.5	12.5	35.0	10.0	3.5	11.0	15.0	390	606UV10B	



TYPE 1FL


Fluorometer cuvette with PTFE cover

• Four windows and base polished



	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
1	3.5	12.5	45.0	1.0	10.0	0.400	1FLG1	1FLUV1	1FLIR1	
2	4.5	12.5	45.0	2.0	10.0	0.700	1FLG2	1FLUV2	1FLIR2	
5	7.5	12.5	45.0	5.0	10.0	1.700	1FLG5	1FLUV5	1FLIR5	
10	12.5	12.5	45.0	10.0	10.0	3.500	1FLG10	1FLUV10	1FLIR10	
20	22.5	12.5	45.0	20.0	10.0	7.000	1FLG20	1FLUV20	1FLIR20	
30	32.5	12.5	45.0	30.0	10.0	10.500	1FLG30	1FLUV30	1FLIR30	
40	42.5	12.5	45.0	40.0	10.0	14.000	1FLG40	1FLUV40	1FLIR40	
50	52.5	12.5	45.0	50.0	10.0	17.500	1FLG50	1FLUV50	1FLIR50	
100	102.5	12.5	45.0	100.0	10.0	35.000	1FLG100	1FLUV100	1FLIR100	

UV 10.0mm

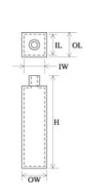

TYPE 11FL

Fluorometer cuvette with quartz cap

• Four windows and base polished

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH WIDTH H		HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
5	7.5	12.5	55.0	5.0	10.0	1.700	11FLUV5	11FLIR5	
10	12.5	12.5	55.0	10.0	10.0	3.500	11FLUV10	11FLIR10	
20	22.5	12.5	55.0	20.0	10.0	7.000	11FLUV20	11FLIR20	
40	42.5	12.5	55.0	40.0	10.0	14.000	11FLUV40	11FLIR40	



TYPE 21FL

IL OL

Fluorometer cuvette with PTFE stopper


• Four windows and base polished

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
5	7.5	12.5	49.0	5.0	10.0	1.700	21FLG5	21FLUV5	21FLIR5	
10	12.5	12.5	49.0	10.0	10.0	3.500	21FLG10	21FLUV10	21FLIR10	
20	22.5	12.5	49.0	20.0	10.0	7.000	21FLG20	21FLUV20	21FLIR20	
40	42.5	12.5	49.0	40.0	10.0	14.000	21FLG40	21FLUV40	21FLIR40	
50	52.5	12.5	49.0	50.0	10.0	17.500	21FLG50	21FLUV50	21FLIR50	





2-2

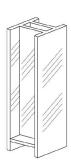


TYPE 52


Dual path length cuvette

• Four windows and base polished

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10 x 1	12.5	12.5	45.0	10.0	1.0	0.350	52UV10X1	52IR10X1	
10 x 2	12.5	12.5	45.0	10.0	2.0	0.700	52UV10X2	52IR10X2	
10 x 3	12.5	12.5	45.0	10.0	3.0	1.050	52UV10X3	52IR10X3	
10 x 4	12.5	12.5	45.0	10.0	4.0	1.400	52UV10X4	52IR10X4	
10 x 5	12.5	12.5	45.0	10.0	5.0	1.750	52UV10X5	52IR10X5	



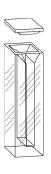

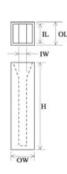


TYPE 9FL

Semi-micro fluorometer cuvette with PTFE cover

- Four windows and base polished
- Base thickness is 3mm

	Outer I	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
10	12.5	12.5	45.0	10.0	4.0	1.400	9FLG10	9FLUV10	9FLIR10	



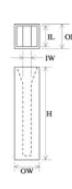

TYPE 18FL

Micro fluorometer cuvette with PTFE cover

- Four windows and base polished
- Base thickness is 3mm


	Outer	Dimensions	(mm)	Inside Dimer		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR
10	12.5	12.5	45.0	10.0	2.0	0.700	18FLG10	18FLUV10	18FLIR10

TYPE 4

Micro fluorometer cuvette with PTFE cover, 5x5mm

• Four windows and base polished

	Outer D	imensio	ns (mm)	Inside Dimen	sions (mm)		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
5x5x33	7.0	7.0	33.0	5.0	5.0	0.640	4G5X5X33	4UV5X5X33	4IR5X5X33	
5x5x45	7.0	7.0	45.0	5.0	5.0	0.800	4G5X5X45	4UV5X5X45	4IR5X5X45	

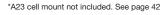



TW TW OW

*A24 cell mount not included. See page 42

TYPE 507

Micro fluorometer cuvette with PTFE cover, 3x3mm

• Four windows and base polished

	Outer D	imensio	ns (mm)	Inside Dimen	sions (mm)	Material & Part Numbers					
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR		
3x3x24	5.0	5.0	24.0	3.0	3.0	0.160	507G3X3X24	507UV3X3X24	507IR3X3X24		
3x3x40	5.0	5.0	40.0	3.0	3.0	0.280	507G3X3X40	507UV3X3X40	507IR3X3X40		



*A24 cell mount not included. See page 42


TYPE 607

Micro fluorometer cuvette with PTFE stopper

• Four windows and base polished

	Outer E	Dimensio	ns (mm)	Inside Dimer	sions (mm)		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH WIDTH		CAPACITY (ml)	UV	IR		
3x3x27	5.0 5.0		27.0	3.0	3.0	0.160	607UV3X3X27	607IR3X3X27		
3x3x43	3 5.0 5.		43.0	3.0	3.0	0.280	607UV3X3X43	607IR3X3X43		



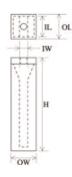


TYPE 29FL

Semi-micro fluorometer cuvette with PTFE stopper

- Four windows polished
- Base thickness is 3mm

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
10	12.5	12.5	49.0	10.0	4.0	1.400	29FLG10	29FLUV10	29FLIR10	



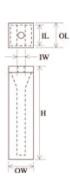

TYPE 30FL

Micro fluorometer cuvette with PTFE stopper

- Four windows polished
- Base thickness is 3mm

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
10	12.5	12.5	49.0	10.0	2.0	0.700	30FLG10	30FLUV10	30FLIR10	






TYPE 701MFL

Sub-micro self masking fluorometer cuvette with PTFE threaded stopper

• Three windows polished

	Outer I	Dimension	s (mm)	Chambe	r Dimensio	ns (mm)			Material & Part Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (µI)	UV
10	12.5	12.5	45.0	10.0	1.0	1.0	8.5	10	701MFLUV10.10A
10	12.5	12.5	45.0	10.0	1.0	1.0	15.0	10	701MFLUV10.10B
10	12.5	12.5	45.0	10.0	1.0	1.0	20.0	10	701MFLUV10.10C
10	12.5	12.5	45.0	10.0	1.0	1.0	8.5	50	701MFLUV10.50A
10	12.5	12.5	45.0	10.0	1.0	1.0	15.0	50	701MFLUV10.50B
10	12.5	12.5	45.0	10.0	1.0	1.0	20.0	50	701MFLUV10.50C
10	12.5	12.5	45.0	10.0	1.0	1.0	8.5	100	701MFLUV10.100A
10	12.5	12.5	45.0	10.0	1.0	1.0	15.0	100	701MFLUV10.100B
10	12.5	12.5	45.0	10.0	1.0	1.0	20.0	100	701MFLUV10.100C
10	12.5	12.5	45.0	10.0	1.0	1.0	8.5	160	701MFLUV10.160A
10	12.5	12.5	45.0	10.0	1.0	1.0	15.0	160	701MFLUV10.160B
10	12.5	12.5	45.0	10.0	1.0	1.0	20.0	160	701MFLUV10.160C

2-4 Screw Cap Fluorometer Cuvettes

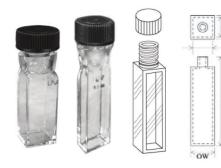

SCREW CAP OD: 16mm

SEPTA: Diameter: 11mm Thickness: 3.5mm

THREADING ID: 8mm THREADING OD: 12mm


TYPE 41FL

Standard fluorometer cuvette with screw cap

- Four windows and base polished
- Each cap contains a PTFE septa seal
- Height includes threading, but not screw cap
- Suitable for anaerobic work and dye laser applications
- Cuvette comes with both open and closed screw caps

	Outer	Dimensions	(mm)	Inside Dimer	sions (mm)		Material & Part Numbers		
L/P (mm)			LENGTH	WIDTH	CAPACITY (ml)	UV			
5	7.5	12.5	55.0	5.0	10.0	1.700	41FLUV5		
10	12.5	12.5	55.0	10.0	10.0	3.500	41FLUV10		
50	52.5	12.5	55.0	50.0	10.0	17.500	41FLUV50		
100	102.5	12.5	55.0	100.0	10.0	35.000	41FLUV100		



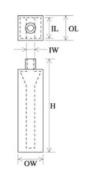

SCREW CAP OD: 16mm SEPTA: Diameter: 11mm Thickness: 3.5mm

THREADING ID: 8mm
THREADING OD: 12mm

TYPE 46FL

Semi-micro fluorometer cuvette with screw cap

- Height includes threading, but not screw cap
- 2 Trought morados timodamig, but not corow dup
- Cuvette comes with both open and closed screw caps
- Suitable for anaerobic work and dye laser applications
- Two windows polished
- Each cap contains a PTFE septa seal
- Base thickness is 3mm

	Outer	Dimensions	(mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV		
10x2	12.5	12.5	55.0	10.0	2.0	0.700	46FLUV10x2		
10x4	12.5	12.5	55.0	10.0	4.0	1.400	46FLUV10x4		








Centrifuge cell

• All surfaces clear fire polished

	Outer Dimens	sions (mm)		Material & Part Numbers	
L/P (mm)	DIAMETER HEIGHT		INNER DIAMETER (mm)	CAPACITY (ml)	UV
7.9	9.9	47.0	7.9	1.700	75UV



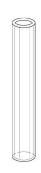


TYPE 76

Scattering cell

- All surfaces clear fire polished
- Bottom is polished

	Outer Dimens	sions (mm)		Material & Part Numbers			
L/P (mm)	" DIAMETER HEIGHT I		INNER DIAMETER (mm)	CAPACITY (ml)	UV		
10	12.0	51.0	10.0	3.100	76UV		



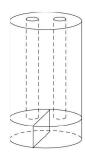

TYPE 528

Cylindrical rotating raman (split) cell

- All surfaces clear fire polished
- Bottom is polished

	Outer Dimens	sions (mm)	Chamber Dime	nsions (mm)	Material & Part Numbers		
L/P (mm)	DIAMETER	HEIGHT	DIAMETER	HEIGHT	CAPACITY (ml)	UV	
19	22.0	35.0	19.0	3.0	0.800	528UV	








TYPE 3FT

Flow-through observation cell, open bottom

- Four windows polished
- Both ends are open

	Outer I	Dimension	s (mm)	Chamber	Dimensio	ns (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	G	UV	IR
10	12.5	12.5	40.0	10.0	10.0	40.0	4.000	3FTG10	3FTUV10	3FTIR10

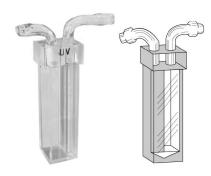


TYPE 45FL

Fluorometer flow cell with detachable inlet/outlet tubes, 2.8ml

- Three windows polished
- Height is without tubes
- Inlet/outlet tubes: OD4 x ID2mm
- Tubes are made of quartz

	Outer I	Dimension	s (mm)	Chambei	Dimensio	ons (mm)		Material & Part Numbers		
L/P (mm)	LENGTH WIDTH HEIGHT		LENGTH WIDTH HEIGHT			CAPACITY UV		IR		
10	12.5	12.5	50.0	10.0	7.0	40.0	2.800	45FLUV10	45FLIR10	



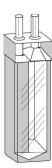

TYPE 58FL

Fluorometer flow-through cell with top tubes

- Three windows polished
- Height is without tubes
- Inlet/outlet tubes: OD4 x ID2 x L10mm Tubes are made of quartz

	Outer I	Dimension	s (mm)	Chamber	r Dimensio	ns (mm)	Material & Part Numbers	
L/P (mm)	LENGTH	WIDTH	HEIGHT			CAPACITY (ml)	UV	
10	12.5	12.5	48.0	10.0	7.0	36.0	2.520	58FLUV10





TYPE 501FL


Fluorometer flow-through cell with top and bottom tubes

- Four windows polished
- Height is without tubes
- Inlet/outlet tubes: OD5 x ID3 x L10mm Tubes are made of quartz

	Outer I	Outer Dimensions (mm) Chamber Dimensions (mm)					Material & Part Numbers	
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	LENGTH WIDTH HEIGHT CAPACITY (ml)		UV	
10	12.5	12.5	65.0	10.0	10.0	33.0	3.300	501FLUV10





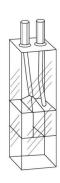

TYPE 59FL

Fluorescence micro flow cell, 440µl

- Four windows and base polished
- Inlet/outlet tubes: OD3 x ID1.5 x L8mm
- Height is without tubes

	Outer I	Dimension	s (mm)	Chamber Dimensions (mm)				Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (ml)	UV	IR
10	12.5	12.5	45.0	10.0	4.0	11.0	15.0	0.440	59FLUV10	59FLIR10




TYPE 74FL

Fluorescence micro flow cell, 30µl

- Four windows polished
- Inlet/outlet tubes: OD2 x ID1
- Height is without tubes

	Outer I	Dimension	s (mm)	Chamber	Dimensio	ns (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	υv	IR
2	3.0	3.0	7.6	2.0	2.0	7.6	0.030	74FLUV2	74FLIR2



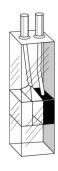

TYPE 514M

Black quartz micro fluorescence flow cell

- Three windows polished
- Inlet/outlet tubes: OD3 x ID1.5 x L10mm
- Height is without tubes
- Suitable for fluorescence and dye laser applications

	Outer I	Dimension	s (mm)	Chamber Dimensions (mm)					Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (ml)	UV	IR	
10	12.5	12.5	45.0	10.0	4.0	11.0	15.0	0.440	514MUV10	514MIR10	






TYPE 608FL

HPLC UV Quartz Flow Cell

- Comes with stainless steel inlet/outlet pipes.
- Dimensions of tubing: Inlet (OD1.6 x ID0.25 x L100mm); Outlet (OD1.6 x ID0.5 x L100mm)
- All cells and tubing tested to 5 bars of pressure (75psi).
- In metal and black quartz housing
- Comes in 3 lightpaths- 1, 1.5 and 3mm
- Available in Z-Dimensions of 8.5 and 15mm
- Three windows optically polished

	Outer Dimensions (mm)			Chamber	r Dimensio	ons (mm)			Material & Part Numbers	
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (μl)	UV	
1	12.5	12.5	45.0	1.0	1.0	11.0	8.5	11	608FLUV1X1X11A	
1	12.5	12.5	45.0	1.0	1.0	11.0	15	11	608FLUV1X1X11B	
1.5	12.5	12.5	45.0	1.5	1.5	11.0	8.5	25	608FLUV1.5X1.5X11A	
1.5	12.5	12.5	45.0	1.5	1.5	11.0	15	25	608FLUV1.5X1.5X11B	
3.0	12.5	12.5	45.0	3.0	3.0	11.0	8.5	100	608FLUV3X3X11A	
3.0	12.5	12.5	45.0	3.0	3.0	11.0	15	100	608FLUV3X3X11B	

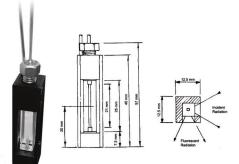

TYPE 8830

μ fluorescence flow cell

• Four windows polished

	Outer I	Dimensions	s (mm)	Chambe	r Dimensio	ns (mm)			Material & Part Numbers	
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (μl)	UV	
1	12.5	12.5	45.0	1.0	1.0	21.0	20.0	20	8830UV1	

This cell has been designed to facilitate the monitoring of HPLC fractions whilst using a spectrofluorometer.


The sample cavity is constructed so that the sample volume is kept to a minimum (20µl nominal). The cell has a large aperture capable of accepting the full exciting and fluorescent radiation. A fluorescence free grade of fused silica has been carefully selected for the construction of the cell which is of a square cross section. The flat sides ensure that any scattering of the incident radiation is kept to a minimum but still allows maximum throughput.

The sample cavity is mounted in an aluminum alloy housing 12.5mm square. Both flow tubes are fitted to the upper face and are 1.5mm OD and 0.25mm ID by 75mm long.

This construction enables the cell to fit any cuvette holder designed for standard fluorescence cuvettes.





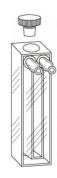

TYPE 54FL

Fluorescence water jacketed cell

- Four windows and base polished
- Tubes to water jacket: OD4 x ID2 x L10mm
- Height is without tubes
- Length is without tubes

2-8

	Outer Dimensions (mm)			Chamber	Dimensio	ns (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR	
10	12.5	12.5	48.0	10.0	4.0	36.0	1.200	54FLUV10	54FLIR10	





TYPE 26FL

Short anaerobic fluorometer cell with glass pouch

- Four windows and base polished
- Overall height is approximately 65mm
- Height is without glass pouch

	Outer Dimensions (mm)			Chambei	r Dimensio	ns (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV		
10	12.5	12.5	45.0	10.0	10.0	30.5	2.500	26FLUV10		



TYPE 28FL

Anaerobic fluorometer cell with glass pouch

- Four windows and base polished
- Overall height is approximately 125mm
- Height is without glass pouch

	Outer I	Dimension	s (mm)	Chambei	r Dimensio	ns (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	G	UV
10	12.5	12.5	84.0	10.0	10.0	30.5	3.500	28FLG10	28FLUV10



TYPE 71FL

Bubble forming anaerobic fluorometer cell with glass pouch

- Four windows and base polished
- Overall height is approximately 102mm
- Height is without glass pouch

	Outer I	Dimension	s (mm)	Chambe	r Dimensio	ns (mm)	Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV	IR	
10	12.5	12.5	76.5	10.0	10.0	45.0	3.500	71FLUV10	71FLIR10	





TYPE 25AFL

Fluorometer Cuvettes

Rectangular tandem fluorometer cell

- Four windows, base, and partition polished
- Height of divider is approximately 37mm

	Outer D	imensior	ns (mm)	Inside Dimen	ısions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH WIDTH HEIGHT			LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10	12.5	12.5	45.0	2 x 4.5	10	2 x 1.300	25AFLUV10	25AFLIR10	



TYPE 56

Rectangular tandem fluorometer cell with PTFE stoppers

• Four windows, base, and partition polished

	Outer E	Dimensior	ns (mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10	12.5	12.5	48.0	2 x 4.5	10	2 x 1.500	56UV10	56IR10	

2-11 Quartz to Glass Grading Fluorometer Cuvettes





TYPE 61FL

Fluorometer cuvette with graded seal tube

- Four windows and base polished
- Quartz to glass graded seal tube: OD8 x ID6mm
- Height is without graded seal tube
- Overall height is approximately 125mm

	Outer D	imensior	ns (mm)	Inside Dimen	nsions (mm)		Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	UV	IR	
10	12.5	12.5	45.0	10	10	3.500	61FLG10	61FLUV10	61FLIR10	




TYPE 62FL

Quartz cuvette with straight bore tube

- All windows polished
- Overall height is approximately 125mm

	Outer D	imensio	ns (mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV		
10	12.5	12.5	140.0	10	10	3.500	62FLUV10		



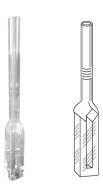

2-11

TYPE 505FL

Semi-micro fluorometer cuvette with graded seal tube

- Four windows and base polished
- Base thickness is 3mm
- Overall height is approximately 125mm
 Quartz to glass graded seal tube: OD8 x ID6mm
- Height is without graded seal tube

	Outer D	imensio	ns (mm)	Inside Dimensions (mm)			Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV		
10	12.5	12.5	45.0	10.0	4.0	1.400	505FLUV10		



TYPE 552

Dual path length cuvette with graded seal tube

- Four windows polished
- Quartz to glass graded seal tube: OD5 x ID3mm
- Height is without graded seal tube
- Overall height is approximately 125mm

	Outer D	imensio	ns (mm)	Inside Dimen	isions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH WIDTH HEIGHT		LENGTH WIDTH		CAPACITY (ml)	UV			
10 x 2	12.5	12.5	45.0	10.0	2.0	0.700	552UV10X2		
10 x 3	12.5	12.5	45.0	10.0	3.0	1.050	552UV10X3		
10 x 4	12.5	12.5	45.0	10.0	4.0	1.400	552UV10X4		
10 x 5	12.5	12.5	45.0	10.0	5.0	1.750	552UV10X5		

2-12 Cryogenic Fluorometer Cuvettes





TYPE 64FL

Cryogenic fluorometer cuvette with PTFE cover

- Four windows and base polished
- Minimum working area of windows is 4mm wide
- All joints are excessively fused to round corners

	Outer D	imensior	ns (mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10	12.5	12.5	45.0	10	10	3.500	64FLUV10	64FLIR10	





TYPE 65FL


Cryogenic fluorometer cuvette with PTFE stopper

- Four windows and base polished
- Minimum working area of windows is 4mm wide
- All joints are excessively fused to round corners

	Outer E	Dimension	ns (mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV		
10	12.5	12.5	49.0	10	10	3.500	65FLUV10		





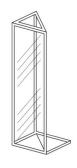


TYPE 81

Triangle cuvette

- Three windows and base polished
- Square base Angle: 45°x45°x90°

	Outer D	imensior	ns (mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10	12.5	12.5	45.0	10	10	1.750	81UV10	81IR10	



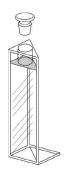

TYPE 82

Triangle cuvette with PTFE stopper


- Three windows and base polished
- Square base Angle: 45°x45°x90°

	Outer D	imensio	ns (mm)	Inside Dimen	nsions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
10	12.5	12.5	49.0	10	10	1.750	82UV10	82IR10	

3-2

Liquid Chromatography and Refractometer Flow Cells



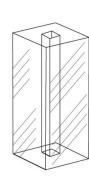

Observation Cell

- Four windows and both ends polished
- Both ends open

TYPE 77

	Outer E	Dimensio	ns (mm)	Inside Dimen	sions (mm)		Material & Part Numbers		
L/P (mm)	LENGTH WIDTH HEIGHT			LENGTH	WIDTH	CAPACITY (ml)			
2	10.5	10.5	27.4	2	2	0.110	77UV2	77IR2	



Cells made from Optical Glass with a wavelength range of 340-2,500nm for use in colorimeters.





TYPE 93

Colorimeter Cell

- Two windows optically polished
- Fits ACS, Data Color, Hunter and several other colorimeters

	Outer I	Dimensio	ns (mm)	Chamber	Dimensi	ons (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	G		
2	6.0	55.0	56.0	2.0	51.0	51.0	5.000	93G2		
5	9.0	55.0	56.0	5.0	51.0	51.0	12.500	93G5		
10	14.0	55.0	56.0	10.0	51.0	51.0	25.000	93G10		
20	24.0	55.0	56.0	20.0	51.0	51.0	50.000	93G20		
30	34.0	55.0	56.0	30.0	51.0	51.0	125.000	93G30		
50	54.0	55.0	56.0	50.0	51.0	51.0	245.00	93G50		









Colorimeter Cell

- Two windows optically polished
- Made for X-Rite colorimeters

	Outer E	Dimensio	ns (mm)	Chamber Dime	ensions (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	
2.5	8.5	28.0	40.0	2.5	24.0	1.800	96G2.5	
5	11.0	28.0	40.0	5.0	24.0	3.600	96G5	
10	16.0	28.0	40.0	10.0	24.0	7.200	96G10	
20	26.0	28.0	40.0	20.0	24.0	14.000	96G20	
40	46.0	28.0	40.0	40.0	24.0	28.000	96G40	
100	106.0	28.0	40.0	100.0	24.0	70.000	96G100	



TYPE 97

Colorimeter Cell

• Two windows optically polished

	Outer D	Dimensio	ns (mm)	Chamber Dime	ensions (mm)	Material & Part Numbers		
L/P (mm)	LENGTH WIDTH HEIGHT		LENGTH	WIDTH	CAPACITY (ml)	G		
40	44.0	20.0	40.0	40.0	16.0	22.400	97G40	




TYPE 529

Small Rectangular Colorimeter

- Two windows optically polished
- Fits ACS, Data Color, Hunter and several other colorimeters

	Outer D	Dimension	ns (mm)	Chamber Dime	ensions (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	G	
5	9.0	24.0	23.20	5.0	20.0	2.32	529G5	



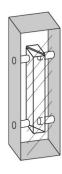

TYPE 510

Refractometer Cell

• Two windows polished

	Outer D	imensio	ns (mm)	Chamber	Dimension	ons (mm)	Material & Part Numbers		
L/P (mm)	LENGTH WIDTH HEIGHT			LENGTH	LENGTH WIDTH HEIGHT			UV	
	5	4.8	12.7	1.6	1.6	8.0	2 x 0.010	510UV	



3-2

TYPE 511

45° Refractometer Flow Cell

• All surfaces and inner chamber windows polished

	Outer D	imensio	ns (mm)	Inside D	imensior	ns (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	LENGTH WIDTH HEIGHT			uv	
	8	10	15	2x2.5	2x1.6	2x7	2x10	511UV	




TYPE 512

15° Refractometer Flow Cell

• All surfaces and inner chamber windows polished

	Outer D	imensio	ns (mm)	Inside D	Dimensior	ns (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (µI)		
	8	10	15	2x2.5	2x1.6	2x7	2x10	512UV	

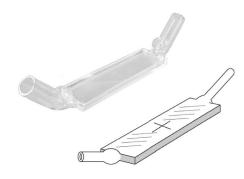

TYPE 78

Electrophoresis cell

- Two windows polished
 - Inlet/outlet tubes: OD7 x ID5 x L20mm
- Height is without tubes

	Outer D	Dimensio	ns (mm)	Inside Dimer	sions (mm)	Material & Part Numbers			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	IR	
1	3.5	12.5	50.0	1.0	10.0	0.475	78UV1	78IR2	
2	4.5	12.5	50.0	2.0	10.0	0.950	78UV2	78IR2	



3-3 **Dye Laser Cells**

TYPE 0005-0023-2

Quanta Ray Dye Laser Flow Through Cell

• Four windows optically polished

	Outer E	Dimensio	ns (mm)	Inside Dimen	sions (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	
6x12	17.5	10.0	16.0	13.5	6.0	12.240	0005-0023-2	





Other Cuvettes

TYPE 0005-0078

Quanta Ray Dye Laser Cell

• Three windows polished

	Outer D	imensio	ns (mm)	Inside Dimen	sions (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	UV	
8	17.5	16.0	10.0	8.0	13.5	10.8	0005-0078	



TYPE 0210

Dye laser flow through cell

- Four windows polished
- Wall thickness is 2mm
- Both ends open
- Originally made for the Quanta Ray[™] dye laser

	Outer E	Dimensior	ns (mm)	Chamber	Dimensi	ons (mm)	Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	uv	
17.0	21.0	17.0	50.8	19.0	15.0	50.8	11.300	0002-0210	





TYPE 2267A

Dye laser cell 20mm

• Four windows and base optically polished

	Outer E	Dimensio	ns (mm)	Inner D	imensior	ns (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	LENGTH WIDTH HEIGHT			uv		
20.0	26.0	26.0	40.0	20.0	20.0	38.0	12.000	2267AUV20		





TYPE 2267B

Dye laser cell 40mm

• Four windows and base optically polished

	Outer D	imensior	ıs (mm)	Inner D	imensior	ıs (mm)		Material & Part Numbers		
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (ml)	UV		
20.0 or 40.0	26.0	46.0	40.0	20.0	40.0	38.0	24.000	2267BUV40		





TYPE 87

Laser absorption cell

- Four windows polished
- Quartz tube: OD6 x ID4 x L20mm
- Height is without tube

	Outer D	Dimensio	ns (mm)	Chamber	Dimensi	ions (mm)		Material & Pa	art Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	CAPACITY (μl)	UV	IR
30.0	32.5	12.5	12.5	30.0	10.0	10.0	3.0	87UV30	87IR30
60.0	62.5	12.5	12.5	60.0	10.0	10.0	6.0	87UV60	87IR60

3-5 Spacers


Spacers are used to reduce volume and light path for standard 10mm l/p cuvettes such as Types 1 & 1FL.

This silica block with a 30mm optical portion is suspended from a flat cover by a fused member which is narrower than the block, to form an overflow chamber when placed in the cuvette.

Spacers cannot be used in a Type 5 cell. Also, please note that due to tolerances, the Type 1 or 1FL you plan to use with the spacer ought to be ordered with the spacer to ensure proper fit.





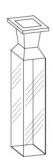

TYPE 527

Precision 2-way insert spacer

• All sides polished

	Oute	er Dimensions (n	nm)	Material & Part Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	uv
1.00 & 5.00	9.00	5.00	46.50	527UV9X5
1.00 & 3.00	9.00	7.00	46.50	527UV9X7
1.00 & 2.00	9.00	8.00	46.50	527UV9X8
1.00 & 0.05	9.00	9.95	46.50	527UV9X9.95
0.50 & 2.00	9.50	8.00	46.50	527UV5X8
0.50 & 1.00	9.50	9.00	46.50	527UV5X9
0.50 & 0.20	9.50	9.80	46.50	527UV9.5X9.8
0.30 & 0.10	9.70	9.90	46.50	527UV9.7X9.9
0.20 & 0.05	9.80	9.95	46.50	527UV9.8X9.95



TYPE 1P

Standard disposable spectrophotometer cuvette

- Two windows clear
- Shipped 100 per box

	Outer E	Dimensio	ns (mm)	Inside Dimer	nsions (mm)			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (μl)	MATERIAL	PART NUMBER
10	12.5	12.5	45.0	10.0	10.0	3.500	UV Plastic	1PUV
10	12.5	12.5	45.0	10.0	10.0	3.500	Acrylic	1PMMA
10	12.5	12.5	45.0	10.0	10.0	3.500	Polystyrene	1PS



TYPE 9P

Semi-micro disposable cuvette

- Two windows clear
- Shipped 100 per box

	Outer Dimensions (mm)			Inside Dimen	mensions (mm)			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (μI)	MATERIAL	PART NUMBER
10	12.5	12.5	45.0	10.0	4.0	1.500	UV Plastic	9PUV
10	12.5	12.5	45.0	10.0	4.0	1.500	Acrylic	9PMMA
10	12.5	12.5	45.0	10.0	4.0	1.500	Polystyrene	9PS



TYPE 21P

Disposable Macro Cuvettes with Stoppers

- Two windows optically polished
- Shipped 100 per box

	Outer D	imensio	ns (mm)	Inside Dimer	sions (mm)		Material & Part Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	POLYSTYRENE
10	12.5	12.5	68.0	10.0	10.0	3.5	21PS



TYPE 509P

Disposable 100mm Tall Macro Cuvettes

- Two windows optically polished
- Shipped 100 per box

	Outer D	imensio	ns (mm)	Inside Dimer	sions (mm)		Material & Part Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (ml)	POLYSTYRENE
10	12.5	12.5	100.0	10.0	10.0	10	509PS





TYPE 700P

Disposable Sub-Micro Cuvettes

- Two windows clear
- Shipped 100 per box

	Outer E	Dimensior	ns (mm)	Inside D	imensio	ns (mm)			Material & Part Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	HEIGHT	ZD (mm)	CAPACITY (µI)	UV PLASTIC
10	12.5	12.5	45.0	10.0	2.0	3.50	8.5	70	700PUV10A
10	12.5	12.5	45.0	10.0	2.0	3.50	15.0	70	700PUV10B



TYPE 704P

Disposable Electroporation Cuvettes

- Two windows clear
- Shipped 100 per box

	Outer E	Dimension	ns (mm)	Inside Dimen	sions (mm)		Material & Part Numbers
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (µl)	POLYCARBONATE
10X1	12.5	12.5	50.0	10.0	1.0	50	704PC10X1
10X2	12.5	12.5	50.0	10.0	2.0	220	704PC10X2
10X4	12.5	12.5	50.0	10.0	4.0	440	704PC10X4



TYPE 1FLP

Fluorometer disposable cuvette

- Four windows clear
- Shipped 100 per box

	Outer D	Dimensior	ns (mm)	Inside Dimen	sions (mm)			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (μI)	MATERIAL	PART NUMBER
10	12.5	12.5	45.0	10.0	10.0	3.500	Acrylic	1FLPMMA
10	12.5	12.5	45.0	10.0	10.0	3.500	Polystyrene	1FLPS



TYPE 9FLP

Fluorometer semi-micro disposable cuvette

- Four windows clear
- Shipped 100 per box

	Outer D	imensio	ns (mm)	Inside Dimer	sions (mm)			
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (µI)	MATERIAL	PART NUMBER
10	12.5	12.5	45.0	10.0	4.0	1.500	Acrylic	9FLPMMA



Disposable Cuvettes, Usable Range

UV PLASTIC						
1PUV	220 ~ 900nm					
9PUV	220 ~ 900nm					

ACRYLIC					
1PMMA	300 ~ 900nm				
1FLPMMA	300nm ~ VIS				
9PMMA	300 ~ 900nm				
9FLPMMA	285 ~ 800nm				

POLYSTYRENE	
1PS	340 ~ 900nm
1FLPS	340 ~ 900nm
9PS	340 ~ 900nm

Compatibility of Disposable Cuvette Materials:

CHEMICAL	POLYSTYRENE	РММА	UV PLASTIC
Acetic acid (96%)	Poor	Poor	Excellent
Acetone	Poor	Poor	Excellent
Acetonitrile	Poor	Poor	Excellent
Benzene	Poor	Poor	Poor
Butanone (methyl ethyl ketone)	Poor	Poor	Excellent
Carbon Tetrachloride	Poor	Poor	Poor
Chloroform	Poor	Poor	Poor
Diethy Ether	Poor	Poor	Poor
Dimethyl Sulfoxide	Poor	Poor	Excellent
Ethanol	Poor	Poor	Excellent
Formaldehye (40%)1	Poor	Poor	Excellent
Heptane	Poor	Poor	Poor
Hexane	Poor	Excellent	Poor
Hydrochloric Acid (32%)	Poor	Poor	Excellent
Isopropanol	Poor	Poor	Excellent
n-Pantane	Poor	Poor	Poor
Oleic Acid	Poor	Poor	Poor
Perchloric Acid (10%)	Poor	Poor	Excellent
Petroleum Ether	Poor	Poor	Poor
Phenol	Poor	Poor	Excellent
Sodium Hydroxide	Excellent	Poor	Excellent
Sulfuric Acid	Poor	Poor	Excellent
Trichloroacetic Acid	Poor	Poor	Excellent
Trichloroethylene	Poor	Poor	Poor
Toluene	Poor	Poor	Poor





PTFE COVERS

DADT#	DESCRIPTION
PART#	DESCRIPTION
A100	PTFE Cover for 100mm cells
	(Type 1-100mm)
A101	PTFE Cover for 1mm cells
A100	(Types 1-1mm, 52-10x1mm)
A102	PTFE Cover for 2mm cells
A10E	(Types 1-2mm, 52-10x2mm)
A105	PTFE Cover for 5mm cells (Types 1 5mm 151 5mm 52 10/5mm)
A110	(Types 1-5mm, 1FL-5mm, 52-10x5mm) PTFE Cover for 10mm cells
ATTU	(Type 1-10mm, 1FL-10mm, 5, 64, 64FL, 508, 1P, 25A, 25AFL)
A120	PTFE Cover for 20mm cells
AIZO	(Type 1-20mm, 1FL-20mm)
A130	PTFE Cover for 30mm cells
ATOU	(Type 1-30mm, 1FL-30mm)
A140	PTFE Cover for 40mm cells
71110	(Type 1-40mm, 1FL-40mm)
A150	PTFE Cover for 100mm cells
	(Type 1-50mm, 1FL-50mm)
A305	PTFE Semi-Micro Cover for 5mm cells
	(Type 9-5mm, 9M-5mm)
A310	PTFE Semi-Micro Cover for 10mm cells
	(Type 9-10mm, 9M-10mm, 9FL-10mm)
A320	PTFE Semi-Micro Cover for 20mm cells
	(Type 9-20mm, 9M-20mm)
A330	PTFE Semi-Micro Cover for 30mm cells
	(Type 9-30mm, 9M-30mm)
A340	PTFE Semi-Micro Cover for 40mm cells
	(Type 9-40mm, 9M-40mm)
A405	PTFE Micro Cover for 5mm cells
	(Type 17-5mm, 18-5mm, 17M-5mm, 18B-5mm, 18BM-5mm, 18M-5mm)
A410	PTFE Micro Cover for 10mm cells
	(Type 17-10mm, 18-10mm, 17M-10mm, 18B-10mm, 18BM-10mm, 18M-10mm)
A420	PTFE Micro Cover for 20mm cells
	(Type 17-20mm, 18-20mm, 17M-20mm, 18B-20mm, 18BM-20mm, 18M-20mm)
A430	PTFE Micro Cover for 30mm cells
	(Type 17-30mm, 18-30mm)
A440	PTFE Micro Cover for 40mm cells
	(Type 17-40mm, 18-40mm)
A450	PTFE Micro Cover for 50mm cells
	(Type 17-50mm, 18-50mm, 18B-50mm)
AT4	PTFE Cover for Type 4
AT507	PTFE Cover for Type 507
AT9B10X1	PTFE Cover for Type 9B-10x1mm
AT9B10X3	PTFE Cover for Type 9B-10x3mm
AT9B10X5	PTFE Cover for Type 9B-10x5mm



PTFE STOPPERS

PART#	DESCRIPTION
B201	PTFE Stopper for 1 and 2mm cells
	(Type 21-1mm, 21-2mm)
B205	PTFE Stopper for 5mm cells
	(Type 21-5mm, 21FL-5mm, 29-5mm, 30-5mm, 82, 29BM 5mm, 29M-5mm, 30BM 5mm,
	30M-5mm, 520M-5mm)
B210	PTFE Stopper for 10mm cells
	(Type 21-10mm, 21-20mm, 21-30mm, 21-40mm, 21-50mm, 21FL-10mm, 21FL-20mm,
	21FL-40mm, 29-10mm, 29-20mm, 29-40mm, 30-10mm, 30-20mm, 30-40mm, 32-10mm,
	35-10mm, 65, 66, 29BM-10mm, 29BM-20mm, 29BM-30mm, 29BM-40mm, 30BM-10mm,
	30BM-20mm, 30M-10mm, 30M-20mm, 515A, 65FL)
B250	PTFE Stopper
	(Type 32-20mm, 34-50mm, 34-100mm, 35-20mm, 35-50mm)
BT50	PTFE Stopper for Type 50
BT53	PTFE Stopper for Type 53
BT54	PTFE Stopper for Type 54 and Type 54FL
BT607	PTFE Stopper for Type 607
D010	PTFE threaded stoppers for sub-micro cells



QUARTZ TOPS

PART#	DESCRIPTION
C010	Quartz Cap (Type 11-10mm, 11-40mm, 11-50mm, 13-10mm, 13-20mm, 13-40mm)
C110	Quartz Stopper for Type 31



SCREW CAPS

PART#	DESCRIPTION
SC1	Closed screw cap (comes with PTFE septa)
SC2	Open screw cap (comes with PTFE septa)
SC3	Mininert Screw Cap Valve
SEPTA	PTFE septa for closed and open screw caps (SC1 and SC2)



SPECIAL TOPS

PART#	DESCRIPTION
A19	M6 Threaded Connectors for Tablet Dissolution Cells (T-601, 602, 605, 606)
A21	Stopcock - Will fit any of the stoppered cells in this catalog.
A22	Funnel - Will fit any of the stoppered cells in this catalog.



These cell mounts will fit into any standard spectrophotometer cell holder. A23 can also be used in fluorometers.

PART#	DESCRIPTION
A20	Cell Mount for rectangular cells less than 5mm light path
	(Type 1-1mm, 1-2mm, 19, 20, 21-1mm, 21-2mm)
A23	Cell Mount suitable for Type 4
A24	Cell Mount suitable for Type 507 and Type 607)
A25	Spacer for 1mm, 2mm and 5mm lightpath cells
A26	Cell mount suitable for Type 42

5-3 Cases & Cell Rack



PART#	DESCRIPTION
CASEA4	Polyfoam case for 45mm high cells - holds 4 cells
CASEB4	Polyfoam case for 49-55mm high cells - holds 4 cells
CASEP2	Plastic (transparent) case to hold 2 cells up to 10mm
CASEP4	Plastic (opaque, black) case to hold 4 cells up to 10mm L/P each, or two 20mm L/P cells,
	or one 50mm L/P cell.
CASE G6	Plastic (opaque, black) case to hold 6 cells up to 10mm L/P each or 3 20mm L/P cells
	or two 50mm L/P cell
P66	12-HOLE CELL RACK: The rack is designed to hold 12 ten millimeter path length cells.
	It is manufactured in polyethylene to ensure that windows cannot be scratched or damaged
	when being placed in the rack.





5-4 Cuvette Washers

PART#	DESCRIPTION
P65S	Single Cell Washer - This model is designed to wash one cell at a time. It comes with a rubber
	stopper for fitting into a size 7 standard vacuum flask. The cell is subjected to a minimum
	amount of handling with less risk of breakage or damage to optical faces. GASKET and
	STOPPER7 included in purchase. (Flask and pump not included.)
P65D	Double Cell Washer - This model is designed to wash two cells at a time and is extremely
	useful for matched pairs of cells, ensuring they are both subjected to the same washing
	treatment. It comes with a rubber stopper for fitting into a size 7 standard vacuum flask.
	The cells are subjected to a minimum amount of handling with less risk of breakage or damage
	to optical faces. GASKET and STOPPER7 included in purchase. (Flask and pump not included.)
GASKET	Small rubber gasket for P65S or P65D.
STOPPER7	Stopper for cell washers to fit size 7 vacuum flask.
Cuvette Washer Full Kit	Full kit for washing cuvettes -contains filtering flask, vinyl tubing, cuvette washer, gasket
	and stopper





CUVETTE CLEANING SOLUTIONS	
MFCSwab	Cleaning swabs for Micro Focus Cell and other Cuvettes
Cuvette Squeegee	Window washer for FireflySci cuvettes

Double-Sided Lens Cloth Double-sided, high quality lens cloth for cleaning optical surfaces on cuvettes

5-5 Mixers/Stirrers





P72





PART#	DESCRIPTION
P67	STIRRING ROD: Manual stirring rod with flared end, manufactured in polyethylene for use in
	10mm cuvettes.
P68	SEMI-MICRO CUVETTE ADD-A-MIXER: Manual mixer fits any semi-micro 10mm pathlength
	cuvette with a 4mm inside with or larger.
P69	CUVETTE ADD-A-MIXER: Manual mixer fits any macro cuvette with 10mm lightpath or larger.
P72	MINI CUVETTE STIRRER: Mini Spectrophotometer Precision Cell Stirrer:
	This PTFE magnetic centrifugal stirrer fits into cuvettes with smaller apertures and volumes
	and provides rapid vertical and horizontal mixing with a minimum of vortexing when placed on
	a magnetic stirring machine. Centrifugal pumping action generated by the cross channels in
	the upper face mixers without aeration. Photometric measurements can proceed uninterrupted
	while reagents are added.
	• Dimensions- 7 x 2mm
	For use with narrow-necked and semi-micro cuvettes
P73	SPECTROPHOTOMETER PRECISION CELL STIRRER: The PTFE magnetic centrifugal stirrer fits
	into the standard 10mm light path spectral cell and provides rapid vertical and horizontal
	mixing with a minimum of vortexing when placed on a magnetic stirring machine. Centrifugal
	pumping action generated by the cross channels in the upper face mixes without aeration.
	Photometric measurements can proceed uninterrupted while reagents are added. 9.5mm
	diameter x 8mm high, the one-piece molded unit will not scratch cuvettes and has no crevices
	to trap dirt. May also be used in test tubes.
P82	ELECTRONIC CELL STIRRER: Stir cell content while measuring absorbance. Fits all standard
	spectrophotometers and fluorometers. The sample is mixed by a tiny magnetic stir bar which
	is placed at the bottom of your cell. A small electronic rotation magnet assembly is placed
	under the cell in the cell compartment. This raises the cell up 5mm which allows you to use
	less sample volume and will not interfere with your unit's operation. Speed of mixer's rotation
	is controlled by an external controller connected by a 36" long paper thin ribbon, allowing your
	instrument to close normally. Controller is 4" x 2.5" x 3" high (including speed controller knob).
	Magnetic stir bar fits any 10x10mm ID cell, and will stir up to 30mm in internal height. The bar

(220-240 VAC Par Number: P82220)

is 8mm high x 9.5mm in diameter. Available for 110-120 or 220-240 VAC wall current.

WHY VALIDATE AND CALIBRATE WITH FIREFLYSCI?

As the world's foremost spectrophotometer reference standard manufacturer, we have a dedicated staff and years of experience to back up our products. Our solid state UV/VIS/NIR reference standards are the most stable on the market and remain in spec long after many other manufacturer's standards go bad. Relying on our proprietary nanodeposition technology, you can rest assured that our solid-state standards will produce steady readings for years to come on your machines.

Other companies keep you bogged down with costly and recurrent recalibrations, which drain time and resources, putting a large drain on your QC department.

The Highest Quality, Most Cost-Effective and Speedy Recalibrations

Let's say that you have some filters from another company already that you'd like to hang onto and recalibrate. No problem. All you need to do is to contact us and let us know. Our average turnaround time is about 1-2 weeks. We will beat anyone's prices and you can proceed knowing that your filters are now backed up by our world class calibration lab.

Advantages at a Glance


- 1. FireflySci is the only company to offer a FREE performance verification service for your current calibration filters.
- 2. Our senior scientists have over 35 years of experience in metrology and quality control.
- 3. Unsurpassed calibration quality.
- 4. Unique products and functionality that no other company offers:
- 5. Combined wavelength/photometric accuracy reference standard (WAV-7 UV/VIS)
- 6. Combined photometric/stray light reference standards (FUV Dual Series)
- 7. Solid-state NIR wavelength reference standard (WAV-8 NIR)
- 8. Continuous photometric accuracy verification in the range of 200-3000nm (SPB Series)
- 9. Metal-on-quartz equivalents work on ANY machine regardless of optics (FUV Series, FUV Dual Series, FNIR-Series, SPB-Series)
- 10. Universal test plate for microplate readers that will work on any company's machines.
- 11. Can use solid state filters for years without having to recalibrate. In over 8 years of active duty, our filters have been proved to never go out of spec.
- 12. Flawless production record: not one return in over 8 years of operation as a result of a manufacturing defect.



Save thousands of dollars a year by using our solid-state photometric accuracy verification calibration standards. Several of our solid state filters are NIST traceable reference standards that never need recalibration due to material aging.

HF -0.5

HF Series: VIS Photometric Accuracy Calibration Standards (400-700nm) (Optical Densities: 0.04 - 4.0au)

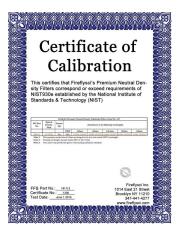

Sample certificate of calibration

HF SERIES STANDARDS

Firefly's HF neutral density solid-state filter series tests photometric accuracy in the VIS range (400-700nm).

PRODUCT SPECS:


- Optical densities available: 0.04au, 0.2au, 0.3au, 0.5au, 0.7au, 1.0au, 1.5au, 2.0au, 2.5au, 3.0au, 4.0au
- Expanded uncertainties of each filter: HF-0.04 (+/- 0.0019au), HF-0.2 (+/- 0.0022au), HF-0.3 (+/- 0.0020au), HF-0.5 (+/- 0.0022au), HF-0.7 (+/- 0.0023au), HF-1.0 (+/- 0.0023au), HF-1.5 (+/- 0.0054au), HF-2.0 (+/- 0.0055au), HF-2.5 (+/- 0.0090au), HF-3.0 (+/- 0.0181au), HF-4.0 (+/- 0.0110au)
- Wavelength range: 400-700nm. Standard validation points included on the certificate of calibration are 440, 465, 546.1, 590, 635nm.
- External dimensions: 12.5x12.5x45mm


BENEFITS OF USING THE HF SERIES PHOTOMETRIC ACCURACY FILTERS:

- Firefly Unique Lifetime Warranty
- Highly accurate and consistent spectrophotometer calibration every time.
- · Never has to be replaced.
- · Never has to be recalibrated if used properly
- Scratch-resistant coating protects optical surfaces
- Can be used in any spectrophotometer (with standard 12.5x12.5x45mm holder)
- Surpasses NIST-930e requirements and tolerances.
- . Full NIST traceability and compatibility.

INCLUDED WITH HF SERIES FILTERS:

- · Instructional manual
- Certificate of calibration with NIST traceability
- Photometric data tables are available for qualified customers. Please contact us for more information.

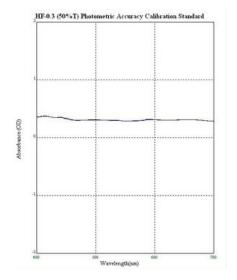

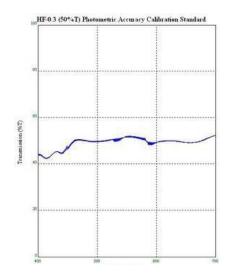

Sample scan charts of HF-Series filter






FUV Series: UV/VIS Photometric Accuracy Calibration Standards (200-700nm) (Optical Densities: 0.04 - 4.0au)

FUV-SERIES

Firefly's FUV neutral density solid-state filter series tests photometric accuracy in the UV and VIS range from 200-700nm.

PRODUCT SPECS:

- Optical densities: 0.04au, 0.2au, 0.3au, 0.5au, 0.7au, 1.0au, 1.5au, 2.0au, 2.5au, 3.0au, 4.0au
- Expanded uncertainties: FUV-0.04 (+/- 0.0019au), FUV-0.2 (+/- 0.0022au), FUV-0.3 (+/- 0.0020au), FUV-0.5 (+/- 0.0022au), FUV-0.7 (+/- 0.0023au), FUV-1.0 (+/- 0.0042au), FUV-1.5 (+/- 0.0054au), FUV-2.0 (+/- 0.0055au), FUV-2.5 (+/- 0.0090au), FUV-3.0 (+/- 0.0181au), FUV -4.0 (+/- 0.0110au).
- Wavelength range: 200-700nm. Standard peaks included on certificate of calibration are 250, 280, 340, 360, 400, 465, 500, 546.1, 590 and 635nm. FireflySci can provide alternate wavelengths in the 200-700nm range upon request.
- External Dimensions: 12.5x12.5x45mm

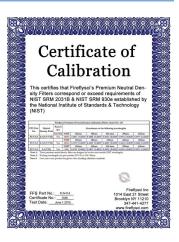

Sample certificate of calibration

BENEFITS:

- Save thousands of dollars on recalibration fees and replacement filters.
- Save time by having standards in your lab where they are needed.
- Save on having to handle fragile and toxic liquid calibration standards.
- Give yourself peace-of-mind knowing that all our standards are 100% NIST Traceable.
- High durability with scratch-resistant coating that protects optical surfaces.
- Can be used in any spectrophotometer in your lab (with standard 12.5 x 12.5 x 45mm holder).



Potassium Dichromate UV/VIS Photometric Accuracy Liquid Standard- NIST 935 (235-430nm)

POTASSIUM DICHROMATE

Potassium Dichromate liquid photometric standards for validating UV/ VIS spectrophotometers.


Potassium Dichromate dissolved in Perchloric Acid (0.001N) is a tried and true method of verifying for absorbance accuracy (photometric) as well as linearity in the UV range. Producing several characteristic peaks throughout the UV range, Potassium Dichromate enables checking at 235nm, 257nm, 313nm, 350nm and 430nm. FireflySci produces its own blend in house, derived from the NIST SRM 935a guidelines and permanently seals the contents in a fire-fused quartz cuvette.

CAPABILITIES:

 \bullet Verify for photometric accuracy and linearity in the UV range at the following wavelengths

235nm 350nm 257nm 313nm

430nm (only available with the 600mg/L concentration)

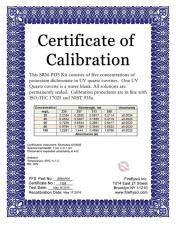

Sample certificate of calibration

INCLUDED WITH POTASSIUM DICHROMATE LIQUID CALIBRATION STANDARD:

- Instructions are included with every Potassium Dichromate liquid filter
- · Certificate of calibration with NIST traceability
- MSDS data for Potassium Dichromate





FNIR Series: NIR Photometric Accuracy Calibration Standards (700-3000nm) (Optical Densities: 0.04 - 4.0au)

FNIR-SERIES


NIST-traceable FNIR series standards for verifying photometric accuracy in the wavelength range of 700-3000nm.

PRODUCT SPECS:

- Optical densities: 0.04au, 0.2au, 0.3au, 0.5au, 0.7au, 1.0au, 1.5au, 2.0au, 2.5au, 3.0au, 4.0au.
- Expanded uncertainties: FNIR-0.04 (+/- 0.0019au), FNIR-0.2 (+/- 0.0022au), FNIR-0.3 (+/- 0.0020au), FNIR-0.5 (+/- 0.0022au), FNIR-0.7 (+/- 0.0023au), FNIR-1.0 (+/- 0.0042au), FNIR-1.5 (+/- 0.0054au), FNIR-2.0 (+/- 0.0055au), FNIR-2.5 (+/- 0.0090au), FNIR-3.0 (+/- 0.0181au), FNIR-4.0 (+/- 0.0110au).
- Wavelength range: 700-3000nm. Certified wavelength points include 700, 1100, 1700, 2700 and 3000nm. Firefly can provide alternative points in the 200-700nm range upon request or full scan data for informational purposes.

• External Dimensions: 12.5x12.5x45mm

• Fully NIST-traceable

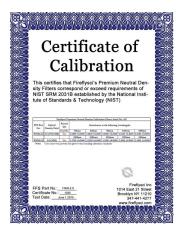

Sample certificate of calibration



Superband Series: Full UV/VIS/NIR
Photometric Accuracy
Calibration Standards
(200-3000nm)
(Optical Densities: 0.04 - 4.0au)

SUPERBAND SERIES


FireflySci is proud to announce the only calibration standard in existence that can verify photometric accuracy over a wavelength range of 200-3000nm.

CAPABILITIES:

 Photometric accuracy in the UV, VIS and NIR range (200-3000nm) at 10 certified wavelengths including: 250, 360, 465, 546.1, 635, 1100, 1700, 2210, 2500, 2850nm.

PRODUCT SPECS:

- Optical densities: 0.04au, 0.2au, 0.3au, 0.5au, 0.7au, 1.0au, 1.5au, 2.0au, 2.5au, 3.0au, 4.0au
- Uncertainties: SPB-0.04 (+/- 0.0019au), SPB-0.2 (+/- 0.0022au), SPB-0.3 (+/- 0.0020au), SPB-0.5 (+/- 0.0022au), SPB-0.7 (+/- 0.0023au), SPB-1.0 (+/- 0.0042au), SPB-1.5 (+/- 0.0054au), SPB-2.0 (+/- 0.0055au), SPB-2.5 (+/- 0.0090au), SPB-3.0 (+/- 0.0181au), SPB-4.0 (+/- 0.0234au).
- Wavelength range: 200-3000nm. FireflySci can provide certified calibration points anywhere in this range. Customers are able to choose points.
- External Dimensions: 12.5x12.5x45mm
- Fully NIST-traceable

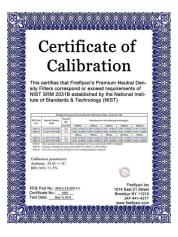

Sample certificate of calibration



NIST Traceable reference standards for verifying the wavelength accuracy of a UV/VIS spectrophotometer. All standards are stringently prepared according to NIST guidelines.

UV/VIS/NIR Calibration Standards

for Spectrophotometers



Liquid Holmium Perchlorate Wavelength Standard- NIST 2034 (240-640nm)

LIQUID HOLMIUM STANDARD


The Firefly Holmium filter is used to calibrate wavelength accuracy in the 240-650 nm range.


CAPABILITIES:

The Firefly Holmium wavelength filter has the following peaks:

- 240.8
 - 451.2
- 277.8
- 485.2
- 287.1 • 361.1
- 536.5 • 640.5
- 416.4

Please note that this is a liquid calibration filter and does not come with the same Lifetime Warranty as Firefly's Solid State Standards. The Firefly liquid holmium filter must be recalibrated every year according to NIST guidelines.

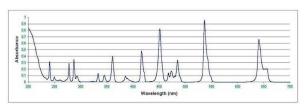

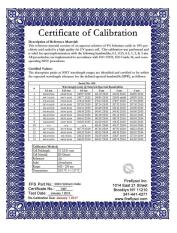

Sample certificate of calibration



WAV-1 UV/VIS Advanced Holmium Oxide Wavelength Accuracy Standard (240-640nm)

WAV-1 UV/VIS


The WAV-1 UV-VIS allows calibration of wavelength accuracy in the ultraviolet (UV) and visible (VIS) range. FireflySci offers the only holmium glass reference standard in the world that has the highly demanded wavelength peak of 241.5nm for verifications relating to DNA measurements and other applications.


CAPABILITIES OF WAV-1 UV/VIS:


We provide the following NIST-traceable peak values for the wavelength accuracy calibration of spectrophotometers:

- 241.5 nm
- 418.5 nm
- 279.3 nm
- 453.4 nm
- 287.6 nm
- 459.9 nm
- 333.8 nm • 360.8 nm
- 536.4 nm • 637.5 nm
- 385.8 nm

Sample absorbance and transmission scans of the WAV-1 UV/VIS:

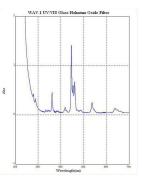

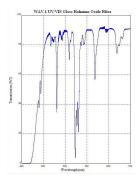

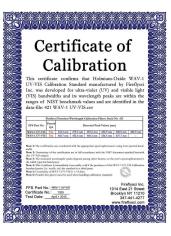

Sample certificate of calibration




WAV-7 UV/VIS Wavelength & UV Photometric Accuracy Standard (329-875nm)


WAV-7

Didymium glass combined wavelength accuracy (329-875nm) and photometric accuracy (270-340nm) filter.


PRODUCT SPECS:

- Wavelength tolerance: +/- 0.2nm
- Wavelength peaks: 329, 472, 512, 681, 875nm
- Photometric tolerance: +/- 0.0024au
- Photometric accuracy verification points: 270, 280, 300, 320, 340nm
- Fully NIST-traceable

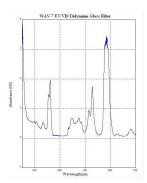

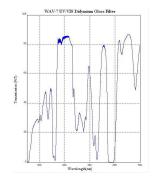

Sample scans of the WAV-7 UV/VIS

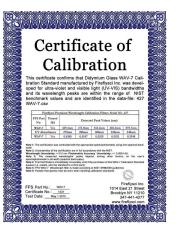

Sample certificate of calibration




WAV-8 NIR Wavelength Accuracy Standard (700-3000nm)

WAV-8


FireflySci's solution for wavelength calibration through the NIR range (700-3000nm).

PRODUCT SPECS:

- Wavelength tolerance: +/- 0.2nm (700-1500nm); +/- 0.5nm (2000-3000nm).
- Approximate wavelength peaks: 707, 747, 803.7, 846, 876.6, 944, 1007, 1518, 1918, 2530, 2719, 3045nm
- Fully NIST-traceable

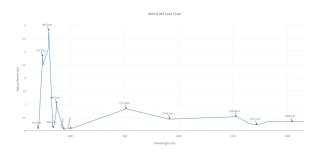

Sample scans of the WAV-8

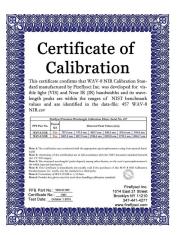

Sample certificate of calibration



Stray Light Calibration

Stray light is defined as light that is detected by the spectrophotometer but does not belong to the given bandwidth of the sample being measured. Often the result of light scattering, diffraction or errors with the instrument, stray light can cause a variety of distortions and errors in analysis.

In order to measure stray light, standards are needed that absorb all light of the wavelength at which the measurement is to be performed and then cut off at a specified wavelength."


UV STRAY LIGHT LIQUID CALIBRATION STANDARDS

UV Stray Light Cut-off Filters ranging from 200-390nm.

CAPABILITIES:

The Firefly Stray Light cut off filters can check for stray light at various wavelength cut-offs in the UV range:

Potassium Chloride: 200nm
Sodium Iodide: 220nm
Potassium Iodide: 260nm
Sodium Nitrite: 390nm

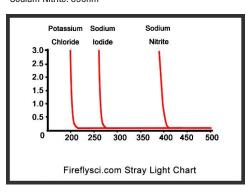



Sample certificate of calibration



UV Stray Light Liquid

Calibration Standards

(200-390nm)


FUV Dual Series: UV/VIS Photometric Accuracy & Stray Light Calibration Standards (200-700nm) (Optical Densities 0.04-4.0au)

FUV DUAL SERIES

Firefly's FUV Dual Series filters allow you to verify photometric accuracy and stray light in the range of 200-700nm. Includes certified values for both the front and back of the filter for determining stray light.

PRODUCT SPECS:

- Optical densities: 0.04au, 0.2au, 0.3au, 0.5au, 0.7au, 1.0au, 1.5au, 2.0au, 2.5au, 3.0au, 4.0au
- Expanded uncertainties: FUV-0.04D (+/- 0.0019au), FUV-0.2D (+/- 0.0022au), FUV-0.3D (+/- 0.0020au), FUV-0.5D (+/- 0.0022au), FUV-0.7D (+/- 0.0023au), FUV-1.0D (+/- 0.0042au), FUV-1.5D (+/- 0.0054au), FUV-2.0D (+/- 0.0055au), FUV-2.5D (+/- 0.0090au), FUV-3.0D (+/- 0.0181au), FUV-4.0D (+/- 0.0110au).
- Wavelength range: 200-700nm. Standard wavelengths included on certificate of calibration are 250, 280, 340, 360, 400, 465, 500, 546.1, 590 and 635nm. FireflySci can provide alternate peaks in the 200-700nm range upon request.
- Includes photometric data for both the front and back side of the filter
- External Dimensions: 12.5x12.5x45mm
- Fully NIST-traceable

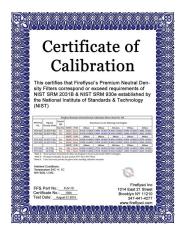

Sample certificate of calibration





NIR Liquid Stray Light Liquid Calibration Standards (800-2000nm)

NIR STRAY LIGHT

A variety of NIST-traceable proprietary blends for checking the stray light of a spectrophotometer in the NIR range from 800-2000nm.

CAPABILITIES:

The Firefly Stray Light cut off filters can check for stray light at various wavelength cut-offs in the NIR range:

- 800nm
- 1000nm
- 1250nm
- 1370nm
- 1550nm
- 1950nm

Spectral Bandwidth Calibration

Spectral Bandwidth (SBW) refers to the physical slit-width of the spectrophotometer's monochromator, which has a direct impact on the resolution of the instruments. Spectral resolution is the ability to resolve two neighboring peaks. Reconciling the resolution is an integral component in getting accurate wavelength and photometric readings on your spectrophotometer. This is a requirement for several pharmacopoeias include European Pharmacopoeia (EP) and United States Pharmacopoeia (USP).

TOLUENE IN HEXANE


PRODUCT SPECS:

Firefly's Toluene in Hexane Spectral Bandwidth Reference Standard is a set containing:

- High-purity solution of 0.02% Toluene in Hexane
- High-purity Hexane blank

SUGGESTED USE:

- The EP recommends comparing the ratio of the peak maxima at 269nm and the minimum at 266nm to determine the spectral resolution of your spectrophotometer.
- Determine the resolution by running the toluene in hexane against the hexane blank between 265-270nm and calculate the peak minima and maxima
- Calculate the ratios of the peak to the trough and compare these results with the reference values.
- Please note that this is a liquid calibration filter and does not come with the same Lifetime Warranty as Firefly's Solid State Standards.
 The Firefly Toluene in Hexane filter must be recalibrated every year according to NIST guidelines.

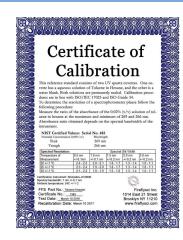

Sample certificate of calibration



Reference kits for NIST, NRC, US Pharmacopoeia and European Pharmacopoeia spectrophotometer validation. All standards are NIST traceable and prepared using NIST guidelines. Custom kits are available upon request.

HF -0.7

NIST 930e VIS Photometric Accuracy Neutral Density Kit (10%T, 20%T, 30%T) (400-700nm)

NIST 930E EQUIVALENT KIT


FireflySci's equivalent of the NIST SRM 930e calibration reference kit for testing photometric accuracy in the VIS range (400-700nm).

THIS KIT INCLUDES THE FOLLOWING FILTERS:

- HF-1 (10%T)
- HF-0.7 (20%T)
- HF-0.5 (30%T)
- · A30 Blank holder
- Case P4

PRODUCT SPECS:

- \bullet Optical densities of the filters: HF-1 (10%T); HF-0.7 (20%T); HF-0.5 (30%T)
- Wavelength range: 400-700nm (this kit can be used for any wavelength within this range)
- Tolerances: HF-0.5 (+/-0.0022au); HF-0.7 (+/- 0.0023au); HF-1.0 (+/- 0.0023au)
- External dimensions of individual filters: 12.5x12.5x45.0mm
- Fully NIST-traceable.

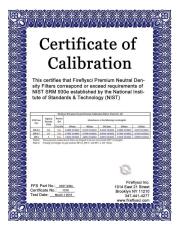

Sample certificate of calibration



NIST 1930 VIS Photometric Accuracy Neutral Density Kit (1%T, 3%T, 50%T) (400-700nm)

NIST 1930 EQUIVALENT KIT


Firefly's equivalent of the NIST SRM 1930 calibration reference kit for testing photometric accuracy in the VIS range (400-700nm).

THIS KIT INCLUDES THE FOLLOWING FILTERS:

- HF-2.0 (1%T)
- HF-1.5 (3%T)
- HF-0.3 (50%T)
- A30 Blank holder
- Case P4

PRODUCT SPECS:

- Wavelength range: 400-700nm (this kit can be used for any wavelength within this range).
- Tolerances: HF-0.3 (+/- 0.0020 au); HF-1.5 (+/-0.0054au); HF-2 (+/-0.0055au)
- External dimensions of individual filters: 12.5x12.5x45.0mm

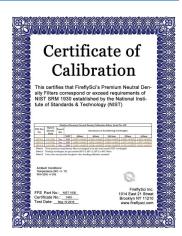

Sample certificate of calibration





NIST 2930 VIS Photometric Neutral Density Kit (0.1%T, 0.3%T, 90%T) (400-700nm)

NIST 2930 EQUIVALENT KIT


Firefly's equivalent of the NIST SRM 2930 calibration reference kit for testing photometric accuracy in the VIS range (400-700nm).

THIS KIT INCLUDES THE FOLLOWING FILTERS:

- HF-3.0 (0.1%T)
- HF-2.5 (0.3%T)
- HF-0.04 (90%T)
- A30 Blank holder
- Case P4

PRODUCT SPECS:

- Wavelength range: 400-700nm (this kit can be used for any wavelength within this range)
- Tolerances: HF-0.04 (+/- 0.0020 au); HF-2.5 (+/- 0.0090au); HF-3 (+/- 0.0181au)
- External dimensions of individual filters: 12.5 x 12.5 x 45.0 mm

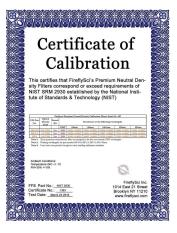

Sample certificate of calibration



Photometric Linearity
Potassium Dichromate Kit
(235-430nm)

PHOTOMETRIC LINEARITY (POTASSIUM DICHROMATE)

Potassium dichromate linearity and photometric accuracy kit for use with many global pharmacopoeia requirements.


PRODUCT SPECS:

- SRM-PD5 Potassium Dichromate Reference Kit contains 5 Concentrations (20mg/L, 40mg/L, 60mg/L, 80mg/L, 100mg/L) and one Perchloric Acid blank.
- SRM-PD6 Kit contains 6 concentrations (20mg/L, 40mg/L, 60mg/L, 80mg/L, 100mg/L, 600mg/L) and one Perchloric Acid Blank. This kit is for use in European Pharmacoepia validations.
- Custom concentrations available on request.

CAPABILITIES:

Verify for photometric accuracy and linearity in the UV range at the following wavelengths

- 235nm
- 257nm
- 313nm
- 350nm
- 430nm (only with the 600mg/L concentration)

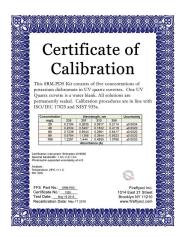

Sample certificate of calibration





NIST 2031 UV/VIS Photometric Accuracy Kit (10%T, 30%T, 90%T) (200-700nm)

NIST 2031 EQUIVALENT KIT

Firefly's upgraded and improved version of the NIST-2031 UV/VIS photometric accuracy kit.

THIS KIT INCLUDES THE FOLLOWING FILTERS:

- FUV-0.04 (90%T)
- FUV-0.5 (30%T)
- FUV-1.0 (10%T)
- A30 Blank Holder
- Case P4

PRODUCT SPECS:

- Optical densities: 0.04au (90%), 0.5au (30%), 1.0au (10%)
- Wavelength Range: 200-700 nm (this kit can be used for any wavelength within this range)
- Tolerances: FUV-0.04 (+/-0.0019 au); FUV-0.5 (+/- 0.0039au); FUV-1 (+/- 0.0042au)
- Wavelength range: 200-700nm. Standard certified calibration points included on the certificate of calibration are 250, 280, 340, 360, 400, 465, 500, 546.1, 590, 635nm. Firefly can provide alternative points in the 200-700nm range upon request or full scan data for informational purposes.
- External Dimensions: 12.5x12.5x45mm



Sample certificate of calibration



FNIR Photometric Accuracy & Linearity Kit (10%T, 20%T, 30%T) (700-3000nm)

FNIR KIT

 $\label{lem:continuous} \text{Kit of three filters for validating photometric accuracy and linearity in the range of } 700-3000 nm$

THIS KIT INCLUDES THE FOLLOWING FILTERS:

- FNIR-1.0 (10%T)
- FNIR-0.7 (20%T)
- FNIR-0.5 (30%T)
- A30 Blank Holder
- Case P4

PRODUCT SPECS:

- Optical densities: 1.0au (10%T), 0.7au (20%T), 0.5au (30%T)
- Wavelength Range: 700-3000nm. Certified wavelength points include 700, 1100, 1700, 2700 and 3000nm. Firefly can provide alternative points in the 700-300nm range upon request or full scan data for informational purposes.
- Tolerances: +/-0.0022 au
- External Dimensions: 12.5x12.5x45mm
- Fully NIST-traceable





Superband UV/VIS/NIR Photometric Accuracy & Linearity Kit (10%T, 30%T, 90%T) (200-3000nm)

SUPERBAND KIT


Kit of three filters and a blank for validating photometric accuracy and linearity in the range of 200-3000nm

THIS KIT INCLUDES THE FOLLOWING FILTERS:

- SPB-0.04 (90%T)
- SPB-0.5 (30%T)
- SPB-1.0 (10%T)
- A30 Blank Holder
- Case P4

PRODUCT SPECS:

- Optical densities: 0.04au (90%T), 0.5au (30%T), 1.0au (10%T)
- Wavelength Range: 200-3000nm. Certified wavelength points include 250, 280, 340, 546.1, 590, 635, 1100, 1700, 2700 and 3000nm. Firefly can provide alternative points in the 200-3000nm range upon request or full scan data for informational purposes.
- Tolerances: +/-0.0022 au
- External Dimensions: 12.5x12.5x45mm

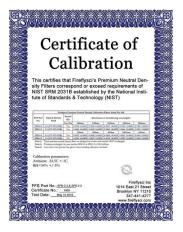

Sample certificate of calibration



VIS Photometric & UV/VIS Wavelength Accuracy Kit

VIS PHOTOMETRIC & UV/VIS WAVELENGTH KIT

Standard calibration kit for checking photometric accuracy in the range of 400-700nm and wavelength accuracy in the range of 241-637nm.

THIS KIT INCLUDES THE FOLLOWING FILTERS:

- HF-0.2 (60%T)
- HF-0.5 (30%T)
- HF-1.0 (10%T)
- WAV-1 UV/VIS Solid State Holmium Oxide Filter
- A30 Blank holder
- Case G6

PRODUCT SPECS:

- Tolerance for photometric accuracy filters: +/- 0.0022 au
- Tolerance for wavelength accuracy filters: +/- 0.2nm
- External dimensions of individual filters: 12.5 x 12.5 x 45.0 mm



European Pharmacopeia Full Compliance Kit

EP KIT

Liquid kit for ensuring that your spectrophotometer is in line with European Pharmacopoeia requirements.

INCLUDED:

- Potassium Dichromate (60mg/L, 600mg/L) with Blank
- Holmium Oxide in Perchloric Acid with Blank
- Potassium Chloride with Pure Water Blank
- Toluene in Hexane with Hexane Blank

PRODUCT SPECS:

- · Check for stray light at wavelength cut-off of 200nm
- Check for spectral bandwidth resolution at 266nm and 269nm
- Check for photometric accuracy at 235nm, 257nm, 313nm, 350nm and 430nm
- Check for wavelength Accuracy at 241.15nm, 287.15nm, 361.5nm and 536.3nm





US Pharmacopoeia Photometric Linearity, Wavelength Accuracy & Stray Light Kit

USP KIT

Liquid kit for ensuring that your UV/VIS spectrophotometer is in full compliance with US Pharmacopeoia standards.

INCLUDED:

- Potassium Dichromate (20mg/L)
- Potassium Dichromate (60mg/L)
- Potassium Dichromate (100mg/L)
- Potassium Dichromate Blank
- Holmium Oxide in Perchloric Acid
- Potassium lodide with Pure Water Blank

CAPABILITIES:

- · Check for stray light at wavelength cut-off of 260nm
- Check for photometric accuracy at 235nm, 257nm, 313nm and 350nm
- Check for wavelength accuracy at 241.15nm, 287.15nm, 361.5nm and 536.3nm



Full UV/VIS Photometric Accuracy, Wavelength & Stray Light Calibration Kit

FULL UV/VIS KIT 4

Broad range solid-state standard kit for checking photometric accuracy, wavelength accuracy and stray light in the UV and VIS range.

THIS CALIBRATION KIT CAN TEST THE FOLLOWING:

- Photometric Accuracy in the UV and VIS range (HF-0.2, HF-0.5, HF-1.0, FUV-0.5)
- Wavelength Accuracy in the UV and VIS range (WAV-1 UV/VIS, WAV-7 UV/VIS)
- Stray Light Performance in the UV and VIS range (FUV-0.5)



Ultimate UV/VIS/NIR
Photometric, Wavelength &
Stray Light Calibration Kit

ULTIMATE UV/VIS/NIR KIT

Broad range solid-state standard kit for checking photometric accuracy, wavelength accuracy and stray light in the UV, VIS AND NIR range.

For use with research-grade UV/VIS/NIR spectrophotometers.

THIS CALIBRATION KIT CAN TEST THE FOLLOWING:

- Photometric Accuracy in the UV, VIS and NIR range (SPB-0.2, SPB-0.5, SPB-1.0)
- Wavelength Accuracy in the UV, VIS an NIR (WAV-1 UV/VIS, WAV-7 UV/VIS, WAV-8 NIR)
- Stray Light Performance in the UV, VIS and NIR range (SPB-1.0D)





Pure Water Fluorescence Reference Standard (250-600nm)



Quinine Sulfate Fluorescence Reference Standard- NIST 936a (375-675nm)

NIST-traceable reference standards for validating fluorescence spectrophotometers

PURE WATER FLUORESCENCE STANDARD

Firefly's ultra-pure water reference standard for validating a fluorescence spectrophotometer

CAPABILITIES:

- Validation of a fluorescence spectrophotometers in the range of 250-600nm.
- Contains two characteristic emission peaks at 397 and 600nm. These are caused by excitations of 350nm and 500nm respectively.

PRODUCT SPECS:

- Firefly Unique Lifetime Warranty
- Highly accurate and consistent calibration every time.
- Can be used in any fluorescence spectrophotometer (with standard 12.5 x 12.5 x 45 mm holder)
- Full NIST traceability and compatibility.

INCLUDED WITH SRM-H20 FLUORESCENCE STANDARD:

- Instruction manual
- · Certificate of calibration with NIST traceability

QUININE SULFATE

Firefly's mixture of Quinine Sulfate in perchloric acid for validating a fluorescence spectrophotometer.

CAPABILITIES:


- Validation of fluorescence spectrophotometers in the range of
- Contains a characteristic emission peak of 450nm by which to calibrate the instrument.

PRODUCT SPECS:

- Firefly Unique Lifetime Warranty
- Comes in a solution of quinine sulfate in 0.105 mol/L of perchloric acid. Comes with an optional perchloric acid blank.
- Highly accurate and consistent calibration every time.
- Can be used in any fluorescence spectrophotometer (with standard 12.5 x 12.5 x 45 mm holder)
- Full NIST traceability and compatibility.

INCLUDED WITH SRM-QS FLUORESCENCE STANDARD:

- · Instruction manual
- · Certificate of calibration with NIST traceability
- . MSDS data for quinine sulfate.

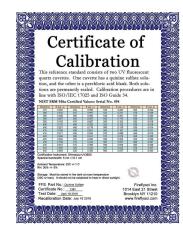


Sample certificate of calibration





Sample certificate of calibration

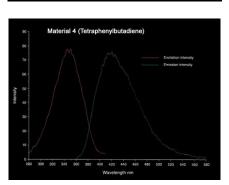

UV/VIS FLUORESCENCE VALIDATION PLATE

UV/VIS fluorescence validation plate for microplate readers.

CAPABILITIES:

• Validation of fluorescence spectrophotometers in the UV/VIS range (290-562nm)


PRODUCT SPECS:


- External dimensions of individual wells: 6.6mm per well
- Total dimensions for 96 well plate: 85x125x12.5mm
- Contains a mix of Anthracene, Compound 610, Napthalene, Ovalene, Terphenyl, Tetraphenylbutadiene, and Rhodamine fluorescent material embedded in PMMA.
- NIST Traceability, ensuring the highest reliability in validation of measurements.


CONTENTS OF FLVP PLATE:


- 1. Wells E1-H1- Mix of **Napthalene** (approximate molar concentration 6X10-5, excitation wavelength 290 nm, emission wavelength 330 nm); corrected excitation 280 nm and emission 330 nm (SBW=5nm) and **Anthracene** (approximate molar concentration 1X10-5, excitation wavelength 360 nm, emission wavelength 402 nm); corrected excitation 360 nm and emission 402 nm (SBW=5nm)
- 2. Wells A7-A10- **p-Terphenyl** (approximate molar concentration 5X10-7, excitation wavelength 295 nm, emission wavelength 338 nm); corrected excitation 295 nm and emission 338 nm (SBW=5nm)
- 3. Wells A1-A4- **Ovalene** (approximate molar concentration 2X10-7, excitation wavelength 342 nm, emission wavelength 482 nm); corrected excitation 342 nm and emission 482 nm (SBW=5nm)
- 4. Wells A12-D12- **Tetraphenylbutadiene** (approximate molar concentration 3X10-7, excitation wavelength 348 nm, emission wavelength 422 nm); corrected excitation 348 nm and emission 422 nm(SBW=5nm)
- 5. Wells H9-H12- **Compound 610** (approximate molar concentration 1X10-6, excitation wavelength 440 nm, emission wavelength 475 nm); corrected excitation 440 nm and emission 480 nm (SBW=5nm)
- 6. Wells H3-H6- **Rhodamine** (approximate molar concentration 2X10-7, excitation wavelength 562 nm, emission wavelength 573 nm); corrected excitation 562 nm and emission 573 nm (SBW=5nm)

Sample scan charts:

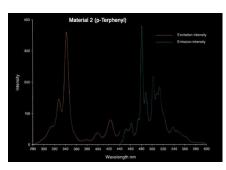

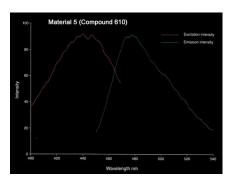

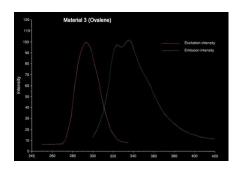

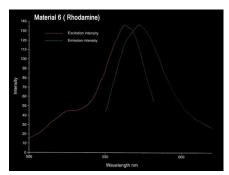

for Spectrophotometers



FireflySci is currently offering a full range of 96 well test plate reference standards including adapters that enable our UV/VIS calibration standards to be used. Test plates cover absorbance, fluorescence and luminescence.

96 WELL PLATE CUVETTE ADAPTER

Cuvette adapter for microplate readers.

CAPABILITIES:

- Enables the user to use any of our cuvette-shaped calibration standards in a 96 well plate form
- You can use the following standards in this adapter: photometric accuracy, wavelength accuracy, stray light, spectral resolution.

PRODUCT SPECS:

- · Made with anodized aluminium for improved performance and resilience.
- Total dimensions for 96 well plate: 85x125x12.5mm
- · Patent-pending design

VIS PHOTOMETRIC & UV/VIS WAVELENGTH VALIDATION PLATE

VIS photometric and wavelength validation plate for microplate readers.

THIS CALIBRATION PLATE INCLUDES THE FOLLOWING WELL STANDARDS (COVERING 16 WELLS IN TOTAL):

- HF-0.2 (60%T)
- HF-0.5 (30%T)
- HF-1.0 (10%T)
- WAV-1 UV/VIS Solid State Holmium Oxide

CAPABILITIES:

- Photometric accuracy in the VIS range (400-700nm)
- Wavelength accuracy in the UV and VIS range (241.5-637.5nm)

PRODUCT SPECS:

- Tolerance for photometric accuracy filters: +/- 0.0022 au
- Tolerance for wavelength accuracy filters: +/- 0.2nm
- · External dimensions of individual wells: 6.6mm per well
- Total dimensions for 96 well plate: 85x125x12.5mm

UV/VIS PHOTOMETRIC & UV/VIS WAVELENGTH ABSORBANCE PLATE

UV/VIS photometric and wavelength validation plate for microplate readers.

THIS CALIBRATION PLATE INCLUDES THE FOLLOWING WELL STANDARDS:

- FUV-0.04 (90%T)
- FUV-0.5 (30%T)
- FUV-1.0 (10%T)
- WAV-1 UV/VIS Solid State Holmium Oxide


CAPABILITIES:

- Photometric accuracy in the UV/VIS range (200-700nm)
- Wavelength accuracy in the UV and VIS range (241.5-637.5nm)

PRODUCT SPECS:

- Tolerance for photometric accuracy filters: +/- 0.0022 au
- Tolerance for wavelength accuracy filters: +/- 0.2nm
- External dimensions of individual wells: 6.6mm per well
- Total dimensions for 96 well plate: 85x125x12.5mm



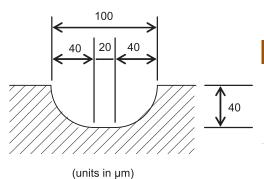

Common Specifications:

- Outer Dimensions: 32 x 65 x t3.8mm
- Material: UV Grade (Synthetic Silica Glass)
- Grooving Method: Wet Etching
- Cross section fluid path (see diagram above)

Most microfluidic chips are custom manufactured. In addition to custom manufacturing, we offer the standard chip configurations as described below. These configurations may offer a starting point upon which modifications can be made. We also offer nanoport assemblies for these chips.

MODEL NAME	FLUID PATH LAYOUT	FLUID PATH LENGTH
YSY-01		28 mm
YLY-01		100 mm
SSS-01	o	55 mm
TST-101		55x20 mm

7-2 Custom Microfluidic Chips


STRUCTURE AND ASSEMBLY

Thermal Fusion: Using special thermal processing, each plate is fused together to create one solid plate. No adhesives are used.

Inlet & Outlet on Sides


Upper plate Fusion

Flow-path plate

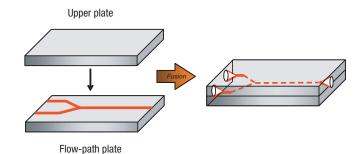



fireflysci

EXAMPLES:

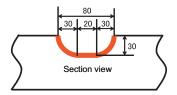

Wet Etching

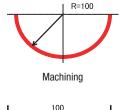

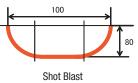


GROOVING METHODS

		Characteristics		
Grooving Method	Width (µm)	Depth (μm)	Surface Roughness (µm)	FEATURES
Wet Etching*	50 ~ 100	5 ~	Ra @ 0.1	High transparency
				Excellent precision in dimensions
				Complicated flow-path layout available
Machining	100 ~ 300	20 ~	Ra @ 0.2	Low cost
				Fair transparency
				Excellent precision in dimensions
				Straight flow-path layout
Shot Blast	100 ~	50 ~	Ra @ 0.7	Low cost
				Complicated flow-path layout available


*WET ETCHING TECHNOLOGY

SHAPE OF ETCHED GROOVE

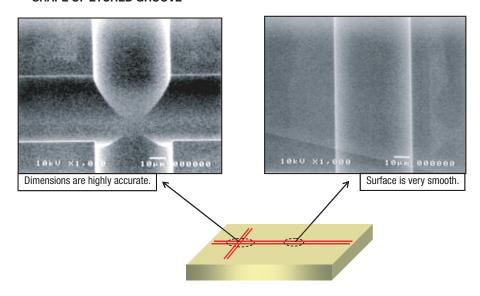

FLOW-PATH CHARACTERISTICS:


Less defects on the edge of the flow path

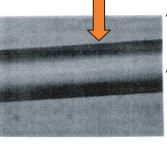

SAMPLE A 60 µm width, 20 µm depth

SAMPLE B 65 µm width, 20 µm depth

No waves on the flow path surface $% \left(1\right) =\left(1\right) \left(1\right$

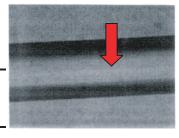



Up to 90 x 90mm is available for the substrate's size.


SUBSTRATE SIZE AND MATERIAL

Our standard material is fused silica (UV Quartz). For details regarding fused silica transmission, absorption, purity, and other properties, please contact us.

?

Please specify groove size (depth and width).



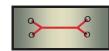

Please specify flow-path layout (examples above).


PROTOTYPE SPECIFICATIONS

In order to manufacture a protoype, please provide the following information:

- 1) Substrate size
- 2) Grooving Method: This will be chosen based on the following:
 - a) Groove size (depth & width)
 - b) Flow-path layout
 - c) Surface roughness
 - d) Specification tolerances
 - e) Required quantity
 - f) Application information (if available)

GENERAL INFORMATION

- FireflySci manufacturing methods guarantee precision down to the smallest micron.
- Production of raw material as well as microchip completion is done within Fireflysci.
- Technical support is available from each level of FireflySci or the life of your product.
- Quality is guaranteed based on ISO 9001 International Standards.

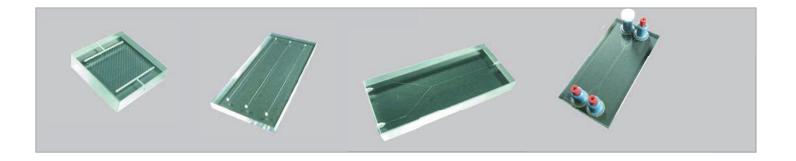





TYPE 526

Standard flow channel cell

• Four windows polished

• Both ends open

	Outer Dimensions (mm)		Inside Dimensions (mm)			Material & Part Numbers	
L/P (mm)	LENGTH	WIDTH	HEIGHT	LENGTH	WIDTH	CAPACITY (μI)	UV
0.25	4.25	4.25	20.00	0.25	0.25	10	526UV0.25

8-2 Custom Flow Channel Cells









Optical Features:

Material: Synthetic Silica Glass Usable Range: 190-2,500nm

Variations:

Single Cone
Double Cone
Polished Cone
Anti-Reflection Coating
Aluminum Coating... and more

Applications:


Particle Counter & Analyzer Blood Cell Counter & Analyzer

Available Size:

Channel Size: 50µm~ Channel Length: 100µm ~ (cone head to cone head)

EXAMPLE OF SOME POSSIBLE MODIFICATIONS

EXAMPLE OF A TYPICAL DESIGN

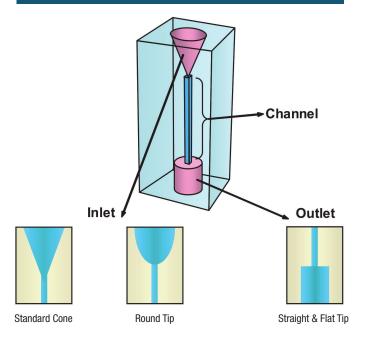

Method: Thermal fusion or optical contact bonding Size: 50 $\mu m \, \sim \,$

Accuracy: \pm 10 µm standard (\pm 5 µm available)

In constructing flow channel cells, four quartz blocks are fused together. Variations using black quartz are also available, as shown above. For a quotation or if you have any questions, please send your drawing or idea to info@fireflysci.com, or fax it to +1-347-554-8048.





Evacuable KBr Dies

KBr Disc Technique: This widely practiced method consists of grinding and mixing a sample with dry KBr or KCl powder, which is then compressed within an evacuated Die to produce a transparent disc, since KBr and KCl are transparent from the ultra-violet well into the infrared. The method may be used for investigation in the ultra-violet and the visible, as well as the infrared region.

Dies suitable for producing discs of 3, 5, 6, 8, 10, 13, and 16mm diameter are available as standard within the size range of the KB series. Non-standard sizes are available on request.

All dies in the KB series are constructed entirely of stainless tool steel, and parts are precision machined. Pellets are optically flat and mirror polished. Tungsten carbide pellets and component parts can be supplied to special order.

Dies are supplied with black anodised sleeves.

Spares for KBr Dies

PART#	DESCRIPTION
KB-30	Set of 1 pellet for KB-3
KB-03	Plunger pellet for KB-3
KB-003	2 sets of 0-rings for KB-3
KB-50	Set of 2 pellets for KB-5
KB-05	Plunger for KB-5
KB-005	2 sets of 0-rings for KB-5
KB-60	Set of 2 pellets for KB-6
KB-06	Plunger for KB-6
KB-006	2 sets of 0-rings for KB-6
KB-80	Set of 2 pellets for KB-8
KB-08	Plunger for KB-8
KB-008	2 sets of 0-rings for KB-8
KB-100	Set of 2 pellets for KB-10
KB-010	Plunger for KB-10
KB-0010	2 sets of 0-rings for KB-10
KB-130	Set of 2 pellets for KB-13
KB-013	Plunger for KB-13
KB-0013	2 sets of 0-rings for KB-13
KB-160	Set of 2 pellets for KB-16
KB-016	Plunger for KB-16
KB-0016	2 sets of 0-rings for KB-16
KB-SE	Spare extractor ring for the above dies
KB-3	KBr Die for 3mm Disc
KB-5	KBr Die for 5mm Disc
KB-6	KBr Die for 6mm Disc
KB-8	KBr Die for 8mm Disc
KB-10	KBr Die for 10mm Disc
KB-13	KBr Die for 13mm Disc
KB-16	KBr Die for 16mm Disc
KB-13ST	KBr Die Starter Kit for 13mm Disc
KB-16ST	KBr Die Starter Kit for 16mm Disc





Evacuable XRF Dies

X-Ray Disc Technique: The BKX series of x-ray dies provides a standard size range producing discs of 20, 25, 32, 35, and 40mm diameter for x-ray fluorescence study of powders and materials. Certain powders which are normally difficult to pelletize without the addition of a binding agent may be successfully pressed within a thin walled plastic ring.

Tungsten carbide plungers and pellets also available.

Spares for X-RF Dies

PART#	DESCRIPTION
KBX-200	Set of 2 pellets for KBX-20
KBX-020	One only plunger for KBX-20
KBX-0020	2 sets of O-rings for KBX-20
KBX-250	Set of 2 pellets for KBX-25
KBX-025	One only plunger for KBX-25
KBX-0025	2 sets of O-rings for KBX-25
KBX-320	Set of 2 pellets for KBX-32
KBX-032	One only plunger for KBX-32
KBX-0032	2 sets of O-rings for KBX-32
KBX-350	Set of 2 pellets for KBX-35
KBX-035	One only plunger for KBX-35
KBX-0035	2 sets of O-rings for KBX-35
KBX-400	Set of 2 pellets for KBX-40
KBX-040	One only plunger for KBX-40
KBX-0040	2 sets of O-rings for KBX-40
KBX-SE	Spare extractor ring for the above dies
KBX-20	Evacuable die for 20mm disc
KBX-25	Evacuable die for 25mm disc
KBX-32	Evacuable die for 32mm disc
KBX-35	Evacuable die for 35mm disc
KBX-40	Evacuable die for 40mm disc



Evacuable 8mm Spanner Die

This KBr Spanner Die is used to produce small numbers of 8mm diameter KBr Discs when a standard die and hydraulic press cannot be justified. The sample is squeezed between the optically polished faces of two 9/16 UNF bolts in the evacuable die body. Moisture is removed from the sample by evacuation. Pressure is then applied by tightening a spanner each end, on the bolts relative to each other. On removal of the bolts, the sample disc remains in the die body. A slide holder is supplied for holding the die body in the FTIR instrument.

PART#	DESCRIPTION
KB/SD-1	Spanner die complete with bolts, seals, & slide holder
KB/SD-2	Spanner die with bolts and seals
KB/SD-3	Spare set of bolts
KB/SD-4	Set of spanners
KB/SD-5	Spare seals
KB/SD-6	Slide holder





KBr Disc Holders with Rectangular Mounting Plate

Other sizes available on request.

PART#	DESCRIPTION
KB-DH3	For 3mm disc
KB-DH5	For 5mm disc
KB-DH10	For 10mm disc
KB-DH13	For 13mm disc
KB-DH16	For 16mm disc
KB-DHM13	Magnetic holder for 13mm disc
KB-FH	Magnetic film holder 20mm aperture



Reflex Evacuable 13mm KBr Minidie (Uni-Die)

Simple, three–piece construction, manufactured from hardened stainless steel. The stepped end pellets are optically polished and there are O-rings on each end for vacuum seal. This can be used with any standard hydraulic press. The pellet is formed in the body and held in the 3" x 2" V-mount for analyzing. After analysis, the sample is washed out and the die is ready for re-use.

PART#	DESCRIPTION	
KB/UD-1	Minidie complete with 3"x2" V-mount	
KB/UD-2	Set of spare stepped pellets (2)	
KB/UD-3	Spare 3"x2" V-mount slide holder	
KB/UD-4	Spare set of 0-rings (4)	



PART#	DESCRIPTION
CHR1/50	Rectangular cuvette holder for cells 1-50mm p/l
CHR1/100	Rectangular cuvette holder for cells 1-100mm p/l
CHC10/50	Cylindrical cuvette holder for cells 10-50mm p/l
CHC10/100	Cylindrical cuvette holder for cells 10-100mm p/l
CHM12.5	Micro cuvette holder for cells 10mm p/l
CHU1/100	Universal cuvette holder for rect & cyl cells 10-100mm

9-3 Gas Cells



100mm Stainless Steel Gas Cell for FTIR

Gas cell is constructed in 316 S.S.

O-rings are used in each end to seal windows with anodized aluminium caps.

Two nupro plug valves are connected via S.S. tubes to body of gas cell.

Window size: 38mm diameter x 6mm, clear aperture 25mm diameter.

Volume is 50cc. Valve outlet connection 1/4" (6.35mm) swagelok.

Two 1/4" (6.35mm) S.S. tubes with rifle ends are supplied for connecting with laboratory tubing.

Slide plate holder with captive screw to lock cell on holder.

Gas cell is supplied with slide plate holder in presentation case.

Windows are not included. Part number: GC-10SS



100mm PTFE Gas Cell for FTIR

Gas cell is constructed in PTFE. O-rings are used in each end to seal windows with anodized aluminium caps.

Two PFA valves are connected to body of gas cell.

Window size: 38mm diameter x 6mm, clear aperture 25mm diameter.

Volume is 50cc. Valve outlet connection 1/4"(6.35mm).

Two 1/4"(6.35mm) PTFE tubes with rifle ends are supplied for connecting with laboratory tubing.

Slide plate holder with captive screw to lock gas cell on holder.

Gas cell is supplied with slide plate holder in presentation case.

Windows are not included. Part number: TGC-10

Product INDEX by Type

50 Cases & Cell Rack	28 T-4	30 T-41FL	38 T-505FL
50 Cell Mounts	07 T- 5	23 T-42	29 T-507
71 Custom Flow Channel Cells	13 T-8	18 T-44	08 T-508
68 Custom Microfluidic Chips	11 T- 9	18 T -45	08 T-509
51 Cuvette Cleaning Solutions	11 T-9B	32 T-45FL	45 T-509P
75 Cuvette Holders for IR & FTIR Instruments	11 T -9BM	17 T-46	41 T-510
50 Cuvette Washers	28 T-9FL	30 T-46FL	42 T-511
63 EP Kit	46 T-9FLP	18 T-48	42 T-512
72 Evacuable KBr Dies	11 T -9M	23 T-49	19 T-513M
73 Evacuable XRF Dies	45 T-9P	20 T-50	33 T-514M
73 Evacuable 8mm Spanner Die	07 T-11	20 T-51H	15 T-520M
64 Full UV/VIS Kit 4	27 T-11FL	20 T-51V	10 T-521
62 FNIR Kit	13 T-16	28 T-52	10 T-522
55 FNIR-Series	13 T-17	21 T-53	09 T-523
54 FUV-Series	13 T-17M	21 T-54	71 T-526
58 FUV Dual Series	14 T-1 8	35 T-54FL	44 T-527
53 HF Series Standards	14 T-18B	37 T-56	31 T-528
74 KBr Disc Holders w/ Rectangular Mounting Plate	14 T-18BM	18 T-58	41 T-529
56 Liquid Holmium Standard	28 T-18FL	32 T-58FL	38 T-552
16 Micro Focus Cell	14 T-18M	33 T-59FL	25 T-601
51 Mixers/Stirrers	22 T-19	19 T-59M	25 T-602
60 NIST 930E Equivalent Kit	23 T-20	24 T-61	26 T-604
60 NIST 1930 Equivalent Kit	08 T-21	37 T-61FL	26 T-605
62 NIST 2031 Equivalent Kit	27 T-21FL	37 T-62FL	26 T-606
61 NIST 2930 Equivalent Kit	45 T-21P	25 T-64	29 T-607
59 NIR Stray Light	22 T-25A	38 T-64FL	34 T-608FL
61 Photometric Linearity	37 T-25AFL	25 T-65	46 T-700P
54 Potassium Dichromate	21 T-26	38 T-65FL	16 T-701M
48 PTFE Covers	36 T-26FL	22 T-71	30 T-701MFL
49 PTFE Stoppers	21 T-28	36 T-71FL	46 T-704P
65 Pure Water Fluorescence Standard	36 T-28FL	19 T-74	34 T-8830
49 Quartz Tops	12 T-29	33 T-74FL	43 T-2267A
65 Quinine Sulfate	12 T-29BM	31 T-75	43 T-2267B
74 Reflex Evacuable 13mm KBr Minidie (Uni-Die)	29 T-29FL	31 T-76	59 Toluene in Hexane
49 Screw Caps	12 T-29M	39 T-77	64 Ultimate UV/VIS/NIR Kit
49 Special Tops	15 T-30	42 T-78	64 USP Kit
68 Standard Microfluidic Chips	15 T-30BM	39 T-81	58 UV Stray Light Liquid Calibration Standards
63 Superband Kit	29 T-30FL	39 T-82	67 UV/VIS Photometric & Wavelength Absorbance Plate
55 Superband Series	15 T-30M	44 T-87	66 UV/VIS Fluorescence Validation Plate
42 T-0005-0023-2	08 T-31	40 T-93	63 VIS Photometric & UV/VIS Wavelength Kit
43 T-0005-0078	09 T-32	40 T-96	67 VIS Photometric & UV/VIS Wavelength Validation Plate
43 T-0210 07 T-1	09 T-34 17 T-34S	41 T-97 20 T-500M	56 WAV-1 UV/VIS 57 WAV-7
27 T-1FL	09 T-35	19 T-500	57 WAV-7
46 T-1FLP	10 T-36	32 T-501FL	67 96 Well Plate Cuvette Adapter
45 T-1P	24 T-37	24 T-504	75 100mm PTFE Gas Cell for FTIR
32 T-3FT	17 T-41	24 T-505	75 100mm Stainless Steel Gas Cell for FTIR
OE I OI I	1 71	∠ 1 000	7 S.I. TOOTHIT CLAIMICOS CLOCK CAS CONTOUT THE


Thank you for your interest in FireflySci, products. We look forward to serving you soon.



About the cuvettes in this catalog...

When you order from FireflySci, Inc., you are ordering from a well established, highly regarded cuvette company. Our cuvettes are 100% guaranteed to meet your quality standards in all regards: Transmission, assembly, matching – everything.

At FireflySci, Inc., quality is number one.

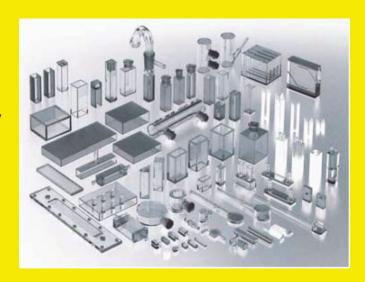

Also inside...

Evacuable dies



FireflySci, Inc.

1014 E. 21st Street, Suite A, Brooklyn, NY 11210 USA

If this catalog did not come with a pricelist, please contact us for one, and take advantage of ordering directly from the manufacturer.

